



Phone: 1-888-347-4537

Fax: 1-724-695-2922

Email: info@clarkehealthcare.com

Web: www.clarkehealthcare.com

Clarke Health Care Products, Inc.

Exclusive Distributors of:

AQUATEC®



ROTH®



AQUATEC®

Ergo VIP

Ergo Dual VIP



**Modular personal care solutions
for complex needs**

ERGONOMIC DESIGN | SAFETY | INDEPENDENCE



Modular personal care solutions for those with more complex seating and positioning needs

Highly customizable


The **Aquatec Ergo VIP** and **Ergo Dual VIP** are modular tilt in space shower chair commodes designed to suit a wide range of individuals and their postural needs. Both models benefit from a -5° to 40° seat tilt angle, whilst the Dual VIP has an additional backrest recline of 0° to 35° for those with more complex seating and positioning requirements.

As with all Ergo models, the VIP and Dual VIP can be tailored to each individual by making some simple adjustments or by adding some easy to fit, off the shelf accessories. Furthermore, their ergonomic design ensures ease of use for the carer whilst providing optimum safety and comfort to the user.

See full range of accessories online at www.clarkehealthcare.com

NEW Contoured, ergonomic seat

- ▶ Based on the imprints people leave behind when seated to ensure a more upright and stable seating position.
- ▶ A 5° seat tilt to lift the knees into the natural squat position for toileting.
- ▶ Transfer supports at the front to allow the user to position themselves into a comfortable position for toileting and showering.

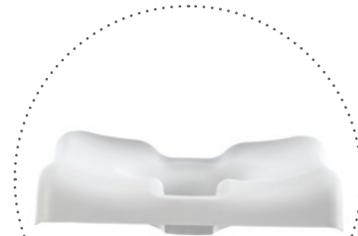

NEW Ergonomic tilt

- ▶ The new pivot and slide mechanism means that during tilt, the chair now moves downwards and forwards which keeps the user's center of gravity close to that of the chair. This makes the tilt operation lighter and more intuitive for the carer, creates a smoother and more comfortable experience for the user and has reduced the overall footprint of the product without affecting its stability.



NEW

FEATURES

-5° to 40°

1

Extremely stable

- ▶ Enhanced stability for additional safety and reassurance

2

Ergonomic tilt

- ▶ Ergonomic tilt with weight shifting abilities

3

Innovative ergonomic seat

- ▶ Contoured, ergonomic seat plate encourages a more upright seating position for enhanced independence and comfort

4

Smaller footprint

- ▶ Optimized wheel base makes it easier to maneuver in small spaces



Other features



▲ Stainless steel corrosion free frame.



▲ Extremely stable even at maximum tilt.



▲ Machine washable tension adjustable backrest.



▲ Armrest height, footrest height and seat height can be easily adjusted.



▲ 5° - 35° Backrest Recline
Ergo Dual VIP Only



Meet the Ergo VIP and Ergo Dual VIP ...

Modular tilt in space models to suit a wide range of individuals and their postural needs.



Ergo VIP

▲ -5°-40° seat tilt angle



Ergo DualVIP

▲ -5°-40° seat tilt angle with an additional backrest recline of 5° - 35°

Technical Data ▼

For more comprehensive information about this product, including the product's user manual, please see your Clarke Health Care rep or Aquatec dealer.



Seat width ▼



Seat depth ▼



Total width ▼



Total depth ▼



Width between armrests ▼

Ergo VIP	18.7"	18.7"	22"	39"	17.75"
Ergo Dual VIP	18.7"	18.7"	22"	39"	17.75"



Tilt angle seat ▼



Tilt angle back ▼



Weight ▼



Load capacity ▼

Ergo VIP	-5° - 40°	-	51.8 lb.	330 lb.
Ergo Dual VIP	-5° - 40°	5° - 35°	62.4 lb.	330 lb.

Frame color ▼



White

Regulatory data ▼




Clarke Healthcare Products
7830 Steubenville Pike, Oakdale, PA 15071
www.clarkehealthcare.com 888-347-4537



AQUATEC®

Ergo Family

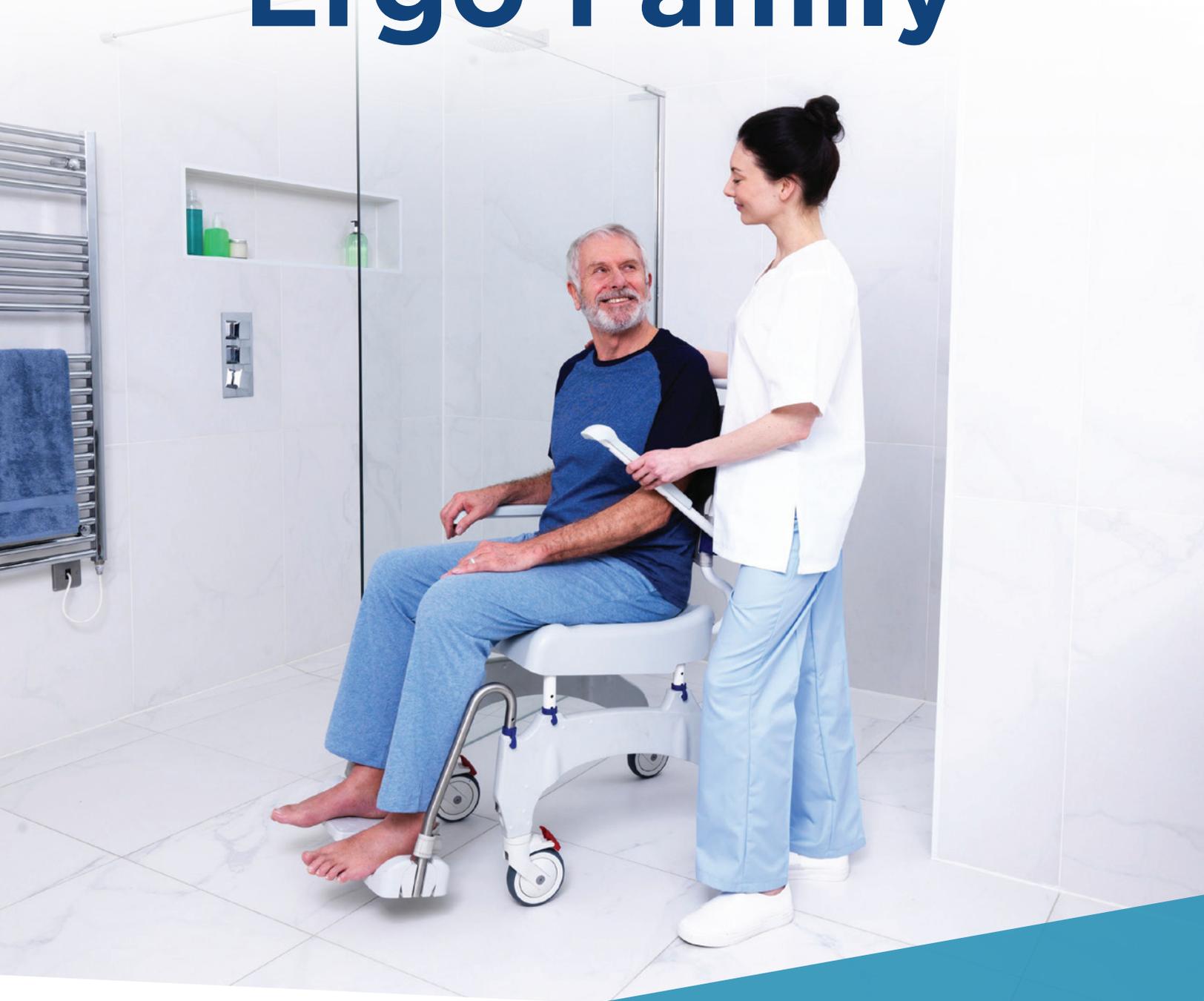

Modular personal care solutions

SAFE | INDEPENDENT | ERGONOMIC



Customer insight-led design to enhance independence, safety and flexibility.

Personal care is highly private activity and is an important part of a person's life and wellbeing. The **Aquatec Ergo Family** of shower chair commodes have been designed with user comfort, independence and dignity in mind. Next generation of the well-established Aquatec shower chairs, the Ergo family comes packed with a range of new features, designed to make personal care activities easier for both the user and carer.

Safety, independence and dignity

A newly contoured seat surface has been developed to encourage a more upright and stable seating position, promoting user independence and overall

comfort. Enhanced stability at the front of the base allows for a smaller footprint and offers user added safety and reassurance when seated in the chair. Furthermore, the optimized wheelbase makes it even easier for a carer to maneuver the chair in small spaces as well as allowing for a better fit over standard toilets for independent toileting.

Modular and flexible design

The modular design of the Ergo family means that the chairs can be adapted to suit a wide range of user needs by making some simple adjustments, or by adding some of our easy to fit, off-the-shelf accessories.

OTHER features include:

- ▶ Stable stainless steel frame
- ▶ Machine washable tension adjustable back with quick-release buckles
- ▶ Easy disassembly for transportation and storage



NEW FEATURES

1 Innovative ergonomic seat

▶ Insight-led design encourages a more upright seating position for enhanced independence and comfort

2 Greater adaptability

▶ With an increased height adjustment range on the footrests and armrests

4 Compact and easy to handle

▶ Optimized wheel base makes it easy to maneuver in small spaces


3 Extremely stable

▶ Enhanced stability at the front for additional safety and reassurance



NEW Contoured, ergonomic seat plate

- ▶ Based on the imprints people leave behind when seated
- ▶ A 5° seat tilt lifts the knees into a natural squat position for toileting
- ▶ Transfer supports at the front to allow the user to position themselves into a comfortable position

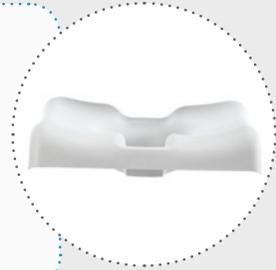

Highly Configurable

Available with a extensive range of positioning, safety, toileting and seating accessories that can be used to support a wide range of conditions and postural needs.

See full range of accessories online at www.clarkehealthcare.com



Positioning



Safety



Toileting



Seating and Backrests

Tilt in space models:

Ergo VIP and Ergo Dual VIP

- ▶ Both have a -5-40° seat tilt angle with a smart weight shift mechanism that ensures optimal stability on a small footprint, whilst the Ergo Dual VIP has an additional backrest recline of 5°-35°



Introducing the **Ergo** Family

From a standard transit model, through to tilt and recline, meet the flexible range of shower chairs that offer personal care solutions for a wide range of physical conditions.



Ergo SPXL ▲

An entry level plus-size model with a maximum user weight of 400 lb.

Ergo XL ▲

An entry level plus-size attendant propelled model with a maximum user weight of 400 lb.

Ergo VIP ▲

Designed for individuals who require slightly more assistance with postural control when seated

Ergo ▲

An entry level attendant propelled model

Ergo Dual VIP ▲

Designed for individuals with complex seating and positioning needs

Ergo SP ▲

An entry level self-propelling model

Technical Data ▼

For more comprehensive information about this product, including the product's user manual, please visit www.clarkehealthcare.com.



Seat width ▼



Seat depth ▼



Total width ▼



Total depth ▼



Width between armrests ▼

	Seat width ▼	Seat depth ▼	Total width ▼	Total depth ▼	Width between armrests ▼
Ergo	18.7"	18.7"	22"	35.3"	17.75"
Ergo XL	18.7"	18.7"	25"	35.3"	20.87"
Ergo SP	18.7"	18.7"	26"	41"	17.75"
Ergo SPXL	18.7"	18.7"	26"	41"	20.87"
Ergo VIP	18.7"	18.7"	22"	39"	17.75"
Ergo Dual VIP	18.7"	18.7"	22"	39"	17.75"



Tilt angle seat ▼



Tilt angle back ▼



Weight ▼



Load capacity ▼

	Tilt angle seat ▼	Tilt angle back ▼	Weight ▼	Load capacity ▼
Ergo	-	-	37.5 lb.	330 lb.
Ergo XL	-	-	39.7 lb.	400 lb.
Ergo SP	-	-	50.7 lb.	330 lb.
Ergo SPXL	-	-	52.9 lb.	400 lb.
Ergo VIP	-5° - 40°	-	51.8 lb.	330 lb.
Ergo Dual VIP	-5° - 40°	5° - 35°	62.4 lb.	330 lb.

Frame color ▼



White

Regulatory data ▼



Clarke Healthcare Products
7830 Steubenville Pike, Oakdale, PA 15071
www.clarkehealthcare.com 888-347-4537



AQUATEC®

Ergo

ACCESSORIES AND ADJUSTMENTS



Accessories

Seating and backrests ▼ ▶



A1603142

- ▶ Standard soft seat overlay. With key shaped cut out.



A1603143

- ▶ Soft seat insert. Can be used with or without the soft seat when commode function is not required.



A1603145

- ▶ Small soft seat overlay. With small, oval opening for smaller and young adults.



A1603152

- ▶ Small soft seat insert. Can be used with the small soft seat when the commode function is not required.



A1603155

- ▶ Splash guard. Clips on easily for use with the standard plastic seat, standard soft seat overlay and small soft seat overlay.



A1531259

- ▶ Ergonomic soft seat overlay, extra soft, textured surface for pelvic stability. 4.6" thick.



A1544469

- ▶ Variable soft seat. Hygiene cut out can be positioned in all directions.



A1603154

- ▶ Universal soft seat overlay. Comfortable closed seat surface.



A1603577

- ▶ XL back. Increases the width between armrests by 3", Included with ErgoXL model, can be retrofitted to standard Ergo model.



A1575299

- ▶ Solid fixed backrest. Offers a stable upright position and is easy to clean.



A1470086

- ▶ Wider XL armrest. Increases distance between armrests by 1.5" each side.



A1558769

- ▶ Soft backrest cushion. For additional stability and comfort. Will also benefit smaller users and young adults. Helps to decrease seat depth.



A1603606, A1604136, A1604016

- ▶ Comfort backrest. For additional comfort, adjustable tension, quick-release buckles. Machine washable.

Positioning ▾



A1535077

- ▶ Side lateral supports. For safe trunk positioning, Width and height adjustable trunk support.



WDSS100P

- ▶ Swing-away side laterals. Adjustable length, swing-away bracket for transfers.



WDACD100

- ▶ Arm trough. Provides additional padded support for arm.



D10167, A1470078

- ▶ Safety Support bar. Additional safety and reassurance, attaches to arm supports. XL size fits ErgoXL.



A1531613

- ▶ Ergonomic headrest holder. Safe and comfortable with various adjustment positions. Pad not included.



A1531608

- ▶ Ergonomic neckrest. Offers exceptional management of the head position.



A1531465

- ▶ Ergonomic headrest. Larger pad for head positioning.



A1470081, A1470082

- ▶ Padded chest belt, seat belt. For safe positioning., Individually adjustable, machine washable.



A1643572, A1643543

- ▶ Calf support with foot support. Padded ELR both height and angle adjustable.



A1643573, A1643574

- ▶ Calf and amputee support. Padded, both height and angle adjustable. For use with the standard plastic seat, standard soft seat overlay and small soft seat overlay.



BB218, BB220

- ▶ Calf strap. Snaps onto footrest tubes, Can prevent feet from kicking behind footrests.



A1470077

- ▶ Arm rest lock. Locks armrest in up or down position.

Toileting ▶



D10230

- ▶ Commode pan with lid. Ergonomic handle allows caretaker to easily remove from pan holder. Scale inside pan.



A1603992, A1601986

- ▶ Front and back anti-tippers. For additional safety in transfers or positioning.

Adjustments

All Ergo models can be adjusted to suit a range of user sizes.

Armrests

- ▶ Ergonomic, swing away armrests can be adjusted to different heights.

Footrests

- ▶ Height adjustable, swing away and detachable footrests
- ▶ No tools required.

Seat Height

- ▶ Seat to ground height can be adjusted between 18.5" - 23.4"
- ▶ No tools required.

Footplates

- ▶ Angle adjustable, flip-up and textured footplates with heel loops
- ▶ No tools required.



Introduction

This useful guide explains how the Aquatec Ergo Family can be customized to suit a wide range of physical conditions by making simple adjustments or by adding some easy to fit, off-the-shelf accessories.

Meet the Ergo Family...



Ergo SPXL ▲

Plus-size self-propelling model

Ergo XL ▲

Plus-size attendant propelled model

Ergo VIP ▲

Tilt-in-space model

Ergo ▲

Attendant propelled model

Ergo Dual VIP ▲

Tilt-in-space model with backrest recline

Ergo SP ▲

Self-propelling model

Compatibility

✓ - Compatible
 X - Not compatible
 S - Included

Most accessories are cross compatible with all Ergo models.

	Product Code	Ergo	Ergo SP	Ergo XL	Ergo SPXL	Ergo VIP	Ergo Dual VIP
Standard soft seat	A1603142	✓	✓	✓	✓	✓	✓
Soft seat insert	A1603143	✓	✓	✓	✓	✓	✓
Small soft seat	A1603145	✓	✓	✓	✓	✓	✓
Small soft seat insert	A1603152	✓	✓	✓	✓	✓	✓
Splash guard	A1603155	✓	✓	✓	✓	✓	✓
Universal Soft seat	A1603154	✓	✓	✓	✓	✓	✓
Ergonomic Special soft seat	A1531259	✓	✓	✓	✓	✓	✓
Variable soft seat	A1544469	✓	✓	✓	✓	✓	✓
Safety bar	D10167	✓	✓	X	X	✓	✓
Safety bar XL	A1470078	X	X	✓	✓	X	X
Chest belt	A1470081	✓	✓	✓	✓	✓	✓
Pelvic/seat belt	A1470082	✓	✓	✓	✓	✓	✓
Armrest lock	A1470077	✓	✓	✓	✓	✓	✓
Arm Trough	WDACD100	✓	✓	✓	✓	✓	✓
XL Back	A1603577	✓	✓	S	S	✓	X
Fixed/solid backrest	A1575299	✓	✓	X	X	✓	✓
Soft backrest cushion	A1558769	✓	✓	✓	✓	✓	✓
Wider XL armrest	A1470086	✓	✓	X	X	✓	✓
24" wheels (set)	A16379	✓	S	✓	S	✓	✓
Calf and Amputee support right	A1556893	✓	✓	✓	✓	✓	✓
Calf and Amputee support left	A1556894	✓	✓	✓	✓	✓	✓
Calf support w/ foot support right ELR	A1546439	✓	✓	✓	✓	✓	✓
Calf support w/ foot support left ELR	A1546705	✓	✓	✓	✓	✓	✓
Calf strap 18"	BB218	✓	✓	✓	✓	✓	✓
Calf strap 20"	BB220	✓	✓	✓	✓	✓	✓
Side lateral	A1535077	✓	✓	✓	✓	✓	✓
Swing-away side lateral	WDSS100P	✓	✓	✓	✓	✓	✓
Ergonomic headrest holder	A1531613	X	X	X	X	✓	✓
Ergonomic headrest	A1531465	X	X	X	X	✓	✓
Ergonomic neckrest	A1531608	X	X	X	X	✓	✓
Commode pan with lid	D10230	✓	✓	✓	✓	✓	✓
Guide rail for commode pan	A1602078	✓	✓	✓	✓	✓	✓
Front and back anti-tipper	A1603992	X	X	X	X	✓	✓
Front anti-tipper	A1601986	✓	✓	✓	✓	X	X
Comfort soft backrest ErgoVIP	A1603606	X	X	X	X	✓	✓
Comfort soft backrest Ergo	A1604136	✓	✓	X	X	X	X
Comfort soft backrest ErgoXL	A1604016	X	X	✓	✓	X	X

Clarke Healthcare Products
 7830 Steubenville Pike, Oakdale, PA 15071
www.clarkehealthcare.com 888-347-4537



AQUATEC® The OCEAN E-VIP



Who is taking care of the caregiver? Aquatec is, with the Ocean E-VIP providing seating and positioning in an adjustable shower commode chair system, designed to alleviate lifting and back strains. Ocean E-VIP consists of a seating system mounted on an electronically height-adjustable wheeled chassis.

For all your hygiene needs, Ocean E-VIP rolls over a commode or into a shower. The seat rises to 41" with a remote control. With the same ease, the remote controls the tilt angle of the seat. The patient is placed at a height or angle best suited for the care necessary. The chair can operate a full day, up to 40 cycles on a single charge. There are plenty of positioning adjustments made with seat, back, head and foot rests. All this movement is done with minimal effort from the caregiver. Ocean E-VIP is suitable for home or institutional use, anywhere a caregiver needs more assistance.

Features

- Adjust seat height by remote
- Tilt-in-space position controlled by remote
- Corrosion-free materials
- Adjustable head support and footrests
- Padded chest belt included
- Directional lock on rear wheel for better handling
- Complete chair can be cleaned in autoclave
- Rechargeable battery provides 40 cycles on a single charge

Components

The frame is constructed from stainless steel, aluminum or plastics to achieve a corrosion free product. The electrical components are sealed, allowing the chair to be cleaned or disinfected by commercial cleaning methods. Recharging can be done on the chair, or when used in shower rooms, the battery case can be removed and recharged off the chair.

Three wheels have foot activated brakes and one rear wheel has a directional lock for better handling when transporting the patient. Padded chest belt is included. The head support can be adjusted horizontally and vertically. Removable, washable back upholstery is also adjustable for best fit.

Warranty

Ocean E-VIP has a limited two year warranty against manufacturer's defects. The battery is excluded.



Sealed electrical components

The OCEAN E-VIP Features



Lowest Height



Raised, Tilted



Seat Raised



Full Height, Tilted



Full Height



Padded Chest Belt



Adjustable Footplate



Directional Caster Lock



Adjustable Head Support

AQUATEC

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

AQUATEC

The OCEAN E-VIP Accessories



Toilet Soft Seat



Action Gel Overlay



Splash Guard



Soft Seat Insert



Safety Bar



Safety Belts



Ergonomic Head Rest Holder



Calf Strap



Calf Support with Foot Support



Calf and Amputee Support



Wide XL Armrests



Collection Pan with Lid



Swing Away Side Lateral



Side Lateral with Height Adjustable Pad



Side Lateral with Angle Adjustable Pad



Arm Trough



Solid Backrest



Soft Backrest Cushion



Head/Neckrest (neckrest shown)



Armrest Locking Device

SPECIFICATIONS

Ocean E-VIP

Overall height, without headrest	41.7" - 63"
Overall width	28.5"
Overall length including footrest	39.4"
Seat height	19.5" - 41"
Seat width	18.75"
Seat depth	17.5"
Width between casters	17.5"
Width between armrests	18"
Width between armrests w/opt. arm width ext.	21"
Battery	24V, 2.9Ah
Weight of chair	86 lb.
Capacity	330 lb.
Product number	A1471208

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

The OCEAN E-VIP

ZITZI™ STARFISH PRO



**NEW
PRODUCT!**

**Combined toilet-
and shower chair**






Zitzi Starfish Pro

A combined shower and toilet chair. Strong construction completely in stainless steel with a lot of adjustment possibilities to provide the best support and position for the user.

Step 1. Hip- & leg position

The combination of a correct set seat part, seat slope, effective hip belt and a good support for the sacrum take away the possibility to shear for the user. This makes it easier to conduct activity and postural control. The adjustment on the seat part is stepless in seat depth, width and abduction/adduction, this allows for a completely individual adaptation of the set part. At the same time Starfish Pro allows the user to grow with the product.

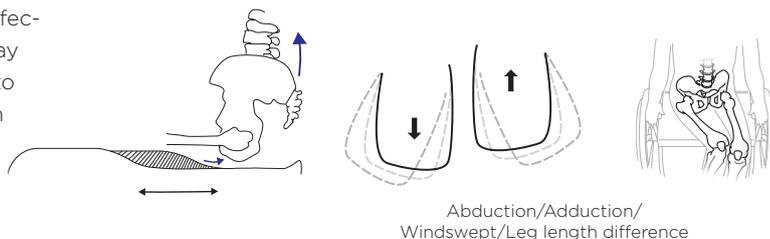

Abduction/Adduction/
Windswept/Leg length difference

Step 2. Placing the hip belt

A correct mounted hip belt positions the pelvic in such a good way so that it prevent shear completely. For Starfish Pro there are many different hip belts available to choose from.




2-points belt
1 mounting per side



4-points belt
2 mountings per side




Anti-rotation belt
2 mountings per side

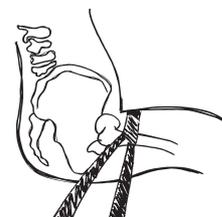

Correct mounted 4 point belt.

Step 3. 3-point positioning

Adjustable trunk support together with hip support create a good possibility to achieve a great positioning for the trunk and this reduces the risk to get in the wrong position and eventually create a malalignment. Good postural control leads to better head control and increases the possibilities for activity and function. Starfish Pro has the possibilities to both correct and keep the users trunk in a good position depending of course on the users possibilities to do so.

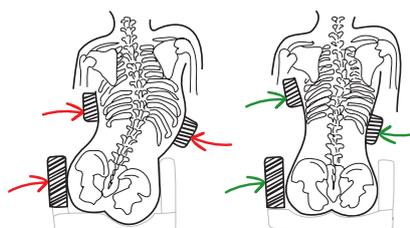


Positioning Harnesses

Together with the hip belt the harness give a good support to be able to position the trunk and make it easier to have a good head control. We provide two different positioning harnesses as a standard option.



Sternum support


Give concentrated support on the sternum.

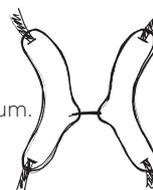

H-strap

Great up stretch for the shoulders and counteract a rotation. Fits great for user with a gavage.

Step 4. Hip angle adjustment

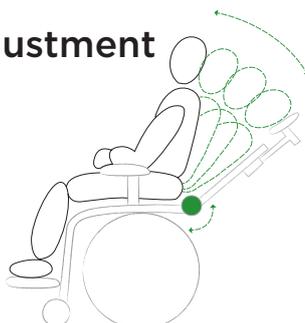

Dynamical supports and materials give freedom of movement, variation and can give a calming effect when You have extension spasticity. The back is angle adjustable which allows the user to get in His/Her perfect hip angle to brake spasticity.

Manual chassi

Headrest

The three membered arm give an optimal adjustability. The headrest is also adjustable in height.

Picture shows:
Zitzi Starfish Pro size 1-2 with adjustable back and driving bow

Positioning

The mounting hardware for the harnesses is adjustable in height and sideways.

Trunk support

Individually adjustable

Height adjusted back part

To increase the postural control

The complete construction is stainless steel.

Manual adjustable height of the seat part.

Footrest

A lot of different footrest to choose from.




Adjustable back

This allows for a better postural control which reduces the risk to get malalignments and increases the possibilities to do activities. In Starfish Pro we are able to both correct and keep the user in a good position depending of course on the users possibilities.



Adjustable tilt and hip angle

Fixed supports is needed to position the user and create a functional seating position. Movement in the materials create a dynamic seat which allows for a certain amount of movement and variation, this creates a calming effect and works great if the user have extension spasticity for example.

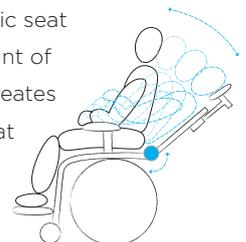

Technical data

USER INFO

User length	46.8-56.5"	54.6-68.2"	66.3-78"
Chassis with seat part	Seat part size 1	Seat part size 2	Seat part size 3
Seat depth	11.7-14.8"	14-18.3"	17.5-21.4"
Seat width	14.8"	17.5"	19.9"
Bredd bendel m foder	14.8"	17.5"	19.9"
	Seat part divided size 1	Seat part divided size 2	Seat part divided size 3
Seat depth	11.7-14.8"	14-18.3"	17.5-21.4"
Seat width	13.6-16(18.7*)"	14.8-17.1(19.9*)"	15.2-17.5(20.2*)"
Width leg art incl cushion	6.4"	7.2"	7.8"
Size hygiene outlet	1.9-4.2(7*)"	1.5-8.9(6.6*)"	1.1-3.5(6.6*)"
Clearance toilet, Man. Chassis	1.6-2"	1.6-2"	1.6-2"
Clearance toilet, Hydraulic/Pwr	1.6-2.5"	1.6-2.5"	1.6-2.5"
Seat height, Man. Chassis	20.6-24.5"	20.6-24.5"	20.6-24.5"
Seat height, Hydraulic/Pwr	20.6-30.4"	20.6-30.4"	20.6-30.4"
Tilt	-5-30°	-5-30°	-5-30°
Hip-/Back angle	-10-35°	-10-35°	-10-30°
Chassis, width	21"	21"	23"
Chassis, length	26.9"	26.9"	30.8"
Wheels	4"	4"	4"
Floor clearance	2.5"	2.5"	2.5"
Ca. weight incl. headrest, footrest, armrest	68.3 lb.	72.7 lb.	77.1 lb.
Max user weight Man. Chassis	265 lb.	265 lb.	330 lb.
Max user weight Hydraulic/Pwr	220 lb.	220 lb.	265 lb.

* Seat parts mounted in the outer tracks

BACK		
Size 1:	Back Height	17.9"
	Back Profile	19.5"
Size 2:	Back Height	21.8"
	Back Profile	23.4"
Size 3:	Back Height	25.7"
	Back Profile	27.3"

LATERAL SUPPORT		
Size 1:	Trunk Width	8.2-11.7"
	Trunk Depth	5.4"
Size 2:	Trunk Width	9.3-13.2"
	Trunk Depth	6.6"
Size 3:	Trunk Width	11.7-15.6"
	Trunk Depth	7.8"
Lateral adjustment +/- 3.3"		

FOOTREST			
	Footrest size 1	Footrest size 2	Footrest size 3
Foot plate (W x D)	13.4 x 8.7"	14.8 x 9.9"	14.8 x 12"
	Footrest divided size 1	Footrest divided size 2	Footrest divided size 3
Foot plate (W x D)	6.2 x 8.5"	7 x 9.9"	7 x 11.9"
Profile Length	11.7"	15.6"	19.5"

High/low adjustable chassi

Headrest

The three membered arm give an optimal adjustability. The headrest is also adjustable in height.



Positioning

The mounting hardware for the harnesses is adjustable in height and sideways.

Trunk support

Individually adjustable

Picture shows:

Zitzi Starfish Pro size 1-2 (divided seat part) with adjustable back and driving bow

Height adjustable back part

To increase the postural control



The complete construction is in stainless steel.



Height adjustable by pump or electric.



Footrest

A lot of different footrest to choose from.




Cover material

Zitzi Starfish Pro series is equipped with a pressure relieving material in the covers. The covers are developed especially for the use in hygiene situations and are both fluid- and dirty fair at the same time it is very pleasant for the skin.

Adjustable leg parts

The leg parts can be adjustable stepless in length, width and in abduction/adduction individually. This is important to be able to get the most optimal seating position for the user at the same time it allows for the user to grow in the product. The seat part is open in the middle this allows for catheterisation.

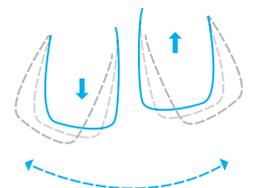

Accessories Starfish Pro



Headrest



Headrest long



Upper mounting brackets



Adjustable back and lateral support



Urine protection



Pot



Table



Calf support



Footrest



Divided footrest swivel adapter



Hip belt



Antirotation belt



Sternum support



H-strap



Bracelet




Ankle holder



Armrest plate size 1




Armrest plate size 2-3

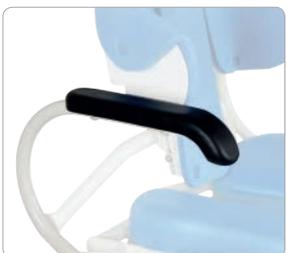

Foldable armrest



Armrest back profile

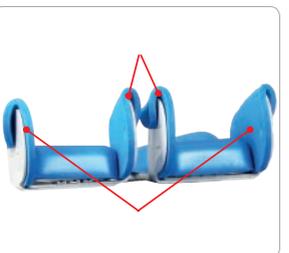

Medial and lateral extensions





Anatomic SITT AB
Box 6137, SE-600 06 Norrköping
Phone +46 11-16 18 00
Email info@anatomicsitt.com
Web www.anatomicsitt.com



Clarke Health Care Products
7830 Steubenville Pike, Oakdale, PA 15071
Phone (724) 695-2122
Email info@clarkehealthcare.com
Web www.clarkehealthcare.com

BORIS...

toilet and shower chair with tilt-in-space positioning

...we achieve more mobility!



Basic model BORIS Size 2

Also available with headrest and foot supports, arm supports and insert pads





The anatomically molded, deep seat pan made of PE plastic with torso and side guides provides particularly high seat stability and tilt-in-space positioning

- Available in five different sizes
- Robust and easy-care sub-frame made of stainless steel
- Height adjustable seat pan for adaptation to the respective constructional situation
- The seat angle is steplessly adjustable via a gas pressure spring from (-) 5° to 30°
- For support in sitting position, the BORIS toilet and shower chair can be combined with the fixation system FLEXI-AKTIV
- We offer head rests and foot supports, insert and arm support pads and a plug-in spray protection as accessories for individual adaptation
- HMV No. 33.40.04.2008



Anatomically molded seat pan with extremely good side guides for stabilization of the upper body. Two-part insert pad for additional sitting comfort



Plug-in arm support for stabilisation and for keeping the upper body upright



Tilt-in-space seat angle is steplessly adjustable via a gas pressure spring from (-) 5° to 30°



Plug-in spray protection available in two sizes



Foldable and angle adjustable foot support made of high-quality ABS material



Height and angle adjustable head rest for support of head control



Practical toilet arrangement

Dimensions inches	Size 0	Size 1	Size 2	Size 3	Size 4
Seat width at front	12	14	16	18	20
Seat width at rear	10	12	14	16	18
Seat depth	12	14	16	18	18.7
Back height	19	20	21.8	23.7	25.3
Seat angle	(-)5 -30°	(-)5 -30°	(-)5 -30°	(-)5 -30°	(-)5 -30°
height over a toilet	18.3-22.2	17.5-21.5	16.3-20.2	16.3-20.2	16.3-20.2
Wide hygiene opening at front	4	4	4	4	4
Wide hygiene opening at back	4.7	4.7	5.3	5.3	5.3
Total width	21.8	23.8	24.5	27.8	28.8
Total length	25.3	25.3	25.3	27.3	27.3
Max. lower leg length	14.8	14.8	18.7	18.7	18.7
Max. load (lb)	132	176	220	220	330

Frame color: Stainless steel



Seatara

A chair for indoor use
Weight 11.6 kg
Self-propelled.

The chair is manufactured using injection technology with advanced plastic substances.

All the metal parts are made from stainless steel (including the bearings).


Hermetic brake locking.

Folding push up handgrips.

Adjustable footrest.

An option for a chamber pot.

An option for a rigid suitcase.

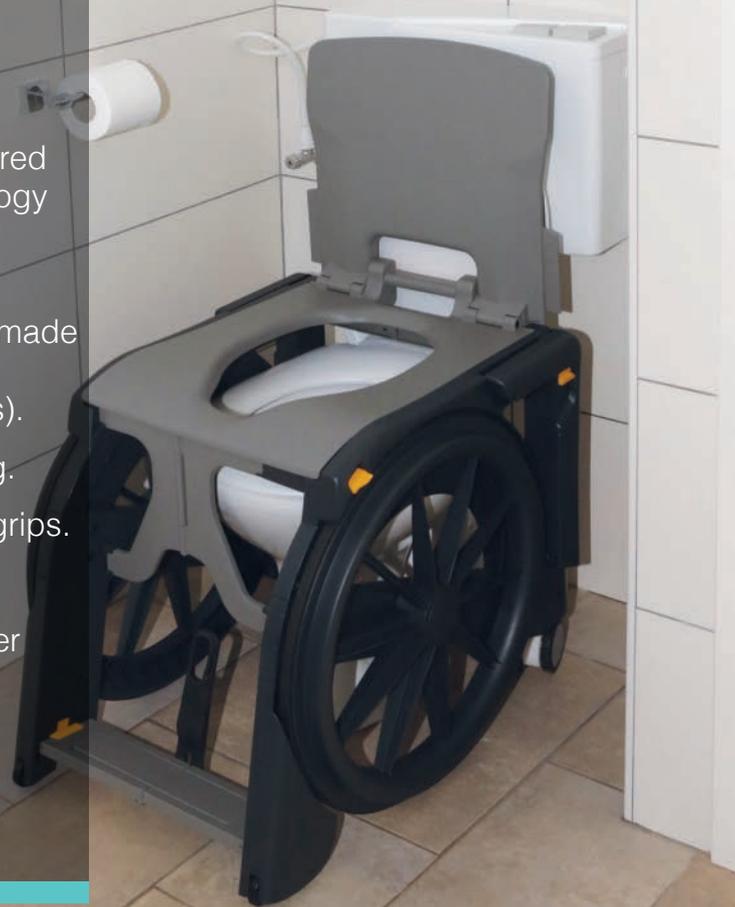

This chair enables safe passage from chair to chair and is convenient for use.



Folding and assembling can be performed by the user while sitting as well.

The small folded dimensions enable the chair's availability for any purpose.



WheelAble

Be Spontaneous

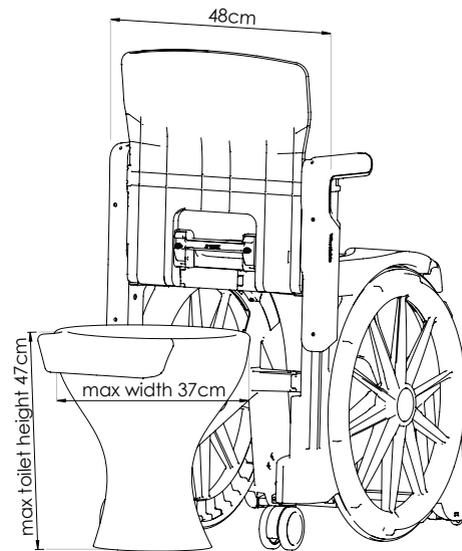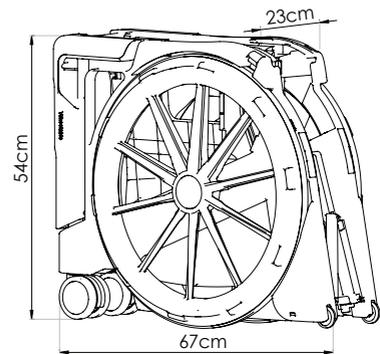

Clarke Health Care Products, Inc
7830 Steubenville Pike
Oakdale, PA 15071
888-347-4537

www.clarkehealthcare.com

Seatarra

WheelAble

Be Spontaneous



A FOLDING CHAIR FOR INACCESSIBLE TOILETS AND SHOWERS

INNOVATION TAKES YOU TO NEW PLACES

Seatar

BathMobile

Designed for you
commode & shower chair



seat turn 180°



adjustable footrest



Weighs 10 kg.
The seat opening can be rotated 180 degrees.
Drive wheels can be assembled.
Long armrests that assist in seating and rising.

Adjustable footrest.

Detachable bedpan rails.

The chair is made with advanced plastic injection technology.

All metal parts are made of stainless steel (including bearings).

features



Seatarra

BathMobile

Designed for you

BathMobile is a breakthrough product in a conservative world.

BathMobile is a modular product that can be detached and assembled without tools, making it easy to move around.

BathMobile allows you the comfort you are used to elsewhere.

BathMobile allows you the freedom to change and adjust; swivel seat, preparation for drive wheels, adjustable product height, detachable bedpan rails.



Distributed by:


Clarke Health Care Products

7830 Steubenville Pike • Oakdale, PA 15071

888-347-4537 • www.clarkehealthcare.com

AQUATEC® The AQUATEC Bathlift

Aquatec shown with Rotary Seat with Transfer Board

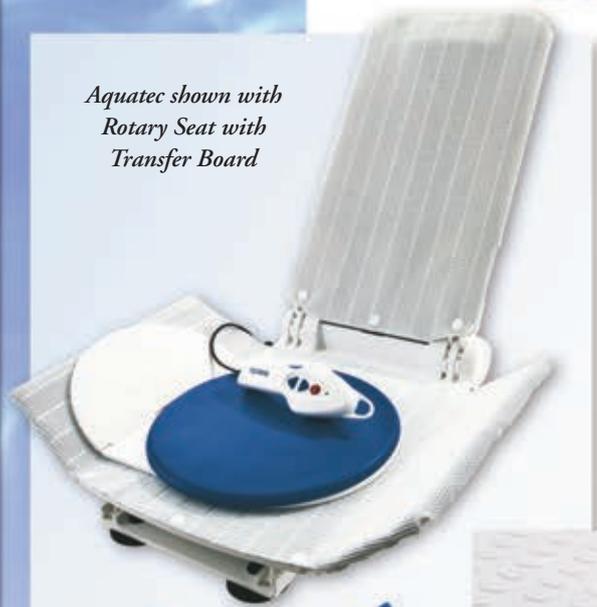

The latest technology for safe and carefree bathing. The **Aquatec** bathlift comes from a family of high quality bath products.

Aquatec is a battery powered bath seat that will lower and raise you again in your tub. Enjoy a long soak and the therapeutic effects of a warm bath. With a reclining back, optional seat widths and back laterals, the **Aquatec** is perfect for those with special bathing needs. The hand remote has easily identified up and down buttons and low power indicator. The lithium-ion battery is in the hand control for convenient recharging. Self release suction cups adhere securely to smooth tub floors and the seat flaps will rest on the top of the tub to allow a stable transfer.

Features

- Backrest with built-in neck support
- Lithium-ion battery with discharge protection, has no memory effect
- Self release suction cups
- Reclining back support up to 40°

Non-slip surface



Components

Aquatec models include a unit to fit standard size tubs, **Aquatec Wide** model for whirlpool style tubs, a back with adjustable side laterals, known as a RSB model and **AquatecXL** model with a 400 pound weight capacity.

All models use a state-of-the-art hand control with integrated battery. The latest lithium-ion manganese-based rechargeable battery pack has total discharge protection and no memory effect. A low battery indicator will light on the hand control when the battery needs to be recharged. A safety feature will not allow the lift to lower when there is not enough power to rise again.

A high backrest with molded neck support can be folded forward by just one hand. The maintenance-free actuator is both powerful and silent. **Aquatec** has a stable and rust-free frame, reinforced lifting scissors and bottom plate.

Aquatec has a non-slip textured surface and comes with removable, washable covers. The bathlift disassembles into easy to carry parts; back, seat and hand control.

AquatecRSB



CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

The AQUATEC Bathlift Accessories



A15461 / A15462
Head support, height adjustable mounts on any back support



A1471210
Seat belt, adjustable cloth strap attaches to seat base


A1471209
Chest belt, used with side laterals, adjustable strap attaches to back support



A13693 / A13694
Sand wedge cushion will stop forward slide associated with high tone or spastic bathers



A14895 / A14901
Abductor block, mounts to lift base

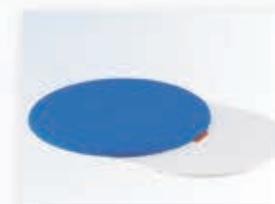

A403002
Rotary seat with transfer board, assists in positioning to center of lift seat. Center disc pivots and slides on base plate



A14586
Carry bag, water resistant nylon with carry handles, holds lift base and back support

AQUATEC

The AQUATEC Bathlift Accessories



A16240
5" Suction pad provides greater area to grab on slightly textured surfaces



A15170 / A14694
Height adapters, adds up to 3 additional inches to base plate, order in 1" increments

Replacement parts including cover mats are also available.

Specifications:


DESCRIPTION	DIMENSIONS	WT. LBS.	THICKNESS
Head Support, Battery Powered Units	12"W x 8"H	1.5	1"
Sand wedge cushion	14"L x 9.5"W x 6"H	16	
Rotary seat with transfer	21"L x 12.5"W	8	.75"
Carry Bag	29"L x 12"W x 14"H	3	
Abductor, A14895	3"Dia. x 4.5"H	.6	
Abductor, A14901	5"Dia. x 5"H	1.6	

AQUATEC

AQUATEC[®] The AQUATEC Bathlift



AquatecRSB Wide

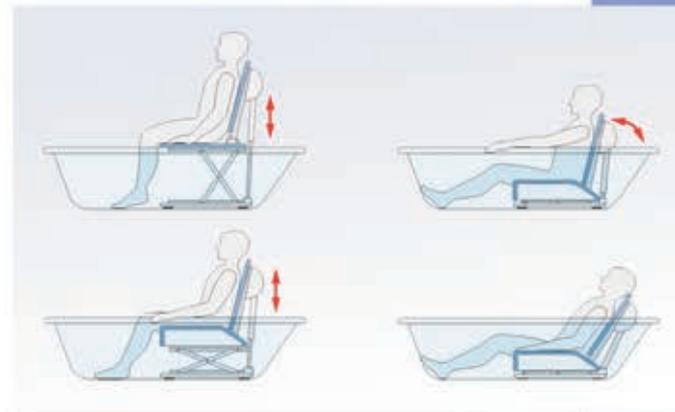

Warranty

Aquatec has a limited three year warranty against manufacturer's defects. Covers, battery and suction cups are excluded.

Specifications:

DESCRIPTION	SPECIFICATION				
	Aquatec	Aquatec Wide	Aquatec RSB	Aquatec RSB Wide	Aquatec XL
Seat width without sideflaps	14.75"	14.75"	14.75"	14.75"	14.75"
Seat width with sideflaps	27.88"	35.5"	27.88"	35.5"	27.88"
Sitting area length	19"	19"	19"	19"	19"
Maximum seat height	17"	17"	17"	17"	17"
Minimum seat height	2.5"	2.5"	2.5"	2.5"	2.5"
Backrest width	14"	14"	27"	27"	14"
Backrest height	26"	26"	26"	26"	26"
Total length (backrest upright)	24"	24"	24"	24"	24"
Total length (backrest reclined)	33"	33"	33"	33"	33"
Total height including backrest	42"	42"	42"	42"	42"
Bottom plate, footprint	12" x 23"	12" x 23"	12" x 23"	12" x 23"	12" x 23"
Cover thickness	.11"	.11"	.11"	.11"	.11"
Total weight without hand control	27 lb	29 lb	31.25 lb	31.25 lb	27 lb
Seat weight	17 lb	19 lb	17 lb	19 lb	17 lb
Backrest weight	10 lb	10 lb	12.5 lb	12.5 lb	10 lb
Hand control weight	1 lb	1 lb	1 lb	1 lb	1 lb
Load capacity	300 lb	300 lb	300 lb	300 lb	400 lb
Backrest reclines to	40°	40°	40°	40°	40°
Part Numbers Blue	A1471565	A1471651	A1471575	A1471582	N/A
Part Numbers White	A1471197	A1471650	A1471202	A1471579	A1471199

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071
Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com



QuattroPlus

Mobile grab handles - the high-load variant.

High carrying load variant: The mobile vacuum grab handle QuattroPlus can be fastened at the wall and across a corner; it serves as horizontal training grab bar and as relocating assistance. The QuattroPlus has been patented in Europe and in the U.S. and can be used also by heavier persons because of its high load carrying capability.

QuattroPlus




Here at Mobeli you can still find them: Experts who after consultation with our technical development department to develop a solution for your very special requirements and who will usually also be able to offer it to you in a short time. Usually our QuattroPlus will be at the heart of a custom design, but there is a multiplicity of combination possibilities with almost all components from the Mobeli system.



A name with substance: The Mobeli system. The high degree of compatibility of the individual grab handle components and the ingeniously simple quick-locking system provide for the possibility to construct whole grab handle railings. We will be pleased to advise you about suitable grab handle combinations for continuous grab handle railings and on the points you should pay particular attention to.



For example you may be looking for the correct Mobeli component system combination for a grab handle railing, e.g. assembled from TÜV-SÜD certified telescopic grab bars with lengths up to 50" and with a nominal carrying load of 150 lbs. You have come to the right place.

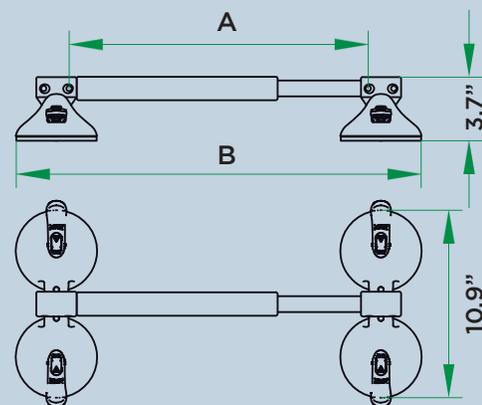

Mobeli building block system



The QuattroPlus as a training grab bar. Because of the integrated cardan joint it is possible to fasten the QuattroPlus at the door frame of a room. A precondition is a pore-tight and smooth surface for attachment, e.g., a metal frame. You only have wooden doors? Inquire about our self-adhesive stainless steel pads.

The QuattroPlus grab handle R1400244SU is designed for door frame internal dimensions from 30.8" to 35.5". By using the diagonal joint adapter these dimensions are extended by 6.5".

As a standard, the integrated cardan joint makes it possible to attach the QuattroPlus grab handle across a corner. If the grip is attached across a corner in the shower, you can pull yourself safely into the shower when entering and hold on to it while taking your shower.



The cardan joints integrated as a standard in QuattroPlus make the QuattroPlus very versatile. The nominal carrying load up to 275 lb. is thereby not impaired.

Integrated cardan joint

QuattroPlus

Mobile grab handles - ingeniously simple and safe.

If the mounting surface is flat, pore-tight, and firm, then the combination of diagonal joint adapter and integrated cardan joint makes nearly any arbitrary attachment possible. In order to withstand the enormous carrying loads that can operate on the grip, the developers at Mobeli have selected a special material combination: aluminium and plastic - durable and versatile.



The diagonal joint adapters specifically developed for the QuattroPlus make it possible to use the QuattroPlus in every direction and up to a nominal carrying load of 176 lb. You can find an explanation of the term "nominal carrying load" on the Mobili brochure on our website.

By using the diagonal joint adapter set, the total length of the grab handle is extended by 6.5".



Thanks to the diagonal joint adapter and its integrated cardan joint, the QuattroPlus can also be mounted diagonally over the bathtub; that's very helpful for getting in and out.



Diagonal joint adapter

MOBELI®	Grip length A	Total length B	Nominal load	Order No.
QuattroPlus	7.8"-12.6"	14"-18.9"	275.6 lb.	R1400241SU
QuattroPlus	12.5"-17.3"	18.7"-24.37"	275.6 lb.	R1400242SU
QuattroPlus	16.3"-21.2"	22.6"-27.49"	242.5 lb.	R1400243SU
QuattroPlus	25"-30"	31.39"-36.27"	245.5 lb.	R1400244SU
Diagonal joint adapter set for QuattroPower			176.4 lb.	R1400557

* See the dimensions of A and B in the sketch on page 2



QuattroPower Support

Power-System - reliable and versatile.

QuattroPower Support

Quick and firm hold without drilling. Just position the QuattroPower supporting system, press down the vacuum tilting lever – ready. You now have a reliable aid for getting up on all smooth and pore-tight surfaces.




The supporting system (QP) QuattroPower Support with its four suction pads distributes the leverage forces acting on the wall. Advantage: It can even be attached on tiled cavity walls (e.g. sandwich plaster board planking), thus eliminating the need for expensive rebuilding to reinforce the walls so that screw-on grab handles might be attached. The QuattroPower Support supporting systems have a nominal carrying load of up to 176 lb.

Very popular near a wall is the combination of a QuattroPower supporting grip with a classical Mobeli grab handle (fixed length or telescopic).



You have a round or mosaic tiled bath tub? Versatile and adaptable, the Mobeli system can offer a solution even in these cases. The QuattroPower supporting grip, mounted upward at the rear tile wall, provides you with safe entry and exit.



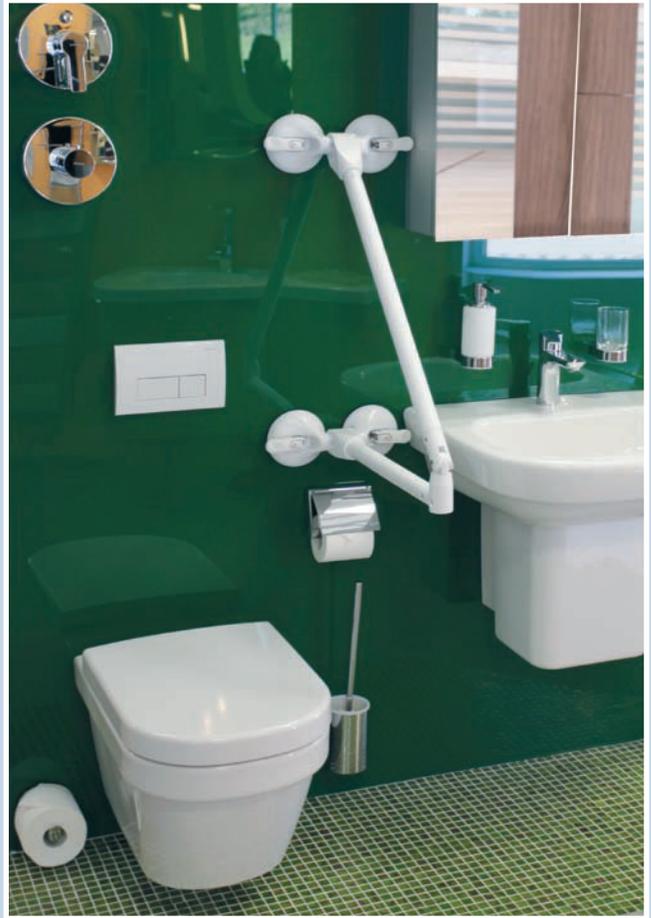

The toilet tissue holder for two rolls of toilet paper can be simply and quickly attached to the QuattroPower.

This handy accessory can also be clipped onto nearly all other mobile grab handles of the Mobeli system.



QuattroPower Support

Depending upon requirements and room conditions, this supporting system can be attached e.g., on both sides of the WC upward and/or downward. An attachment upward offers additional possibilities for support or pulling up. Our QP supporting systems offer a multitude of attachment variants and permit use with partially completed new tasks as, e.g., in physiotherapy.


An additional grasp horn for the QP Support facilitates getting up and can also be used as a toilet tissue holder.



All grips of the QuattroPower supporting system are equipped with swivelling double joints and built-in safety brakes that prevent an inadvertent lifting up and yet enable folding up to the wall.

Well-suited for travelling. You can quickly and simply disassemble the QP. It takes up little space in your luggage and is quickly operational again when needed.

Space-saving minimal solution: You can rebuild the QP to a simple grab handle with four suction pads and attach it directly at the wall or across a corner.

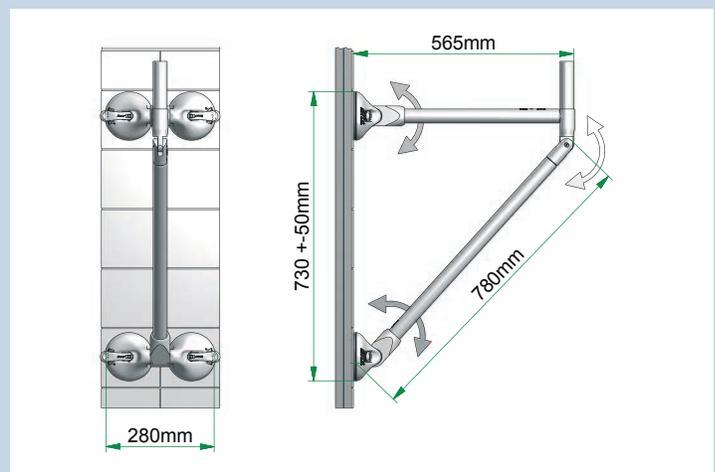

MOBELI®	Grip length	Nominal load	Order No.
QuattroPower Support	22"	176 lb.	R1400265SU
QuattroPower Support with grip horn	22"	176 lb.	R1400266SU
Dual toilet paper holder			R1400650



The Mobeli supporting system QuattroPower Tub can be attached inside and outside of the bathtub. It is able to endure high pull and vertical loads and even lateral loads of up to 110 lb. of nominal carrying load.

Good to know: The suction pads hold even under water.

If you do not need the supporting system any longer you can fold it up. A joint stopper holds the lifted QuattroPower Tub in the vertical position.



QuattroPower Tub

Power system - reliable and versatile.

Optimized for the bathtub: The Mobeli supporting system QuattroPower Tub is excellently suited for getting into and out of the bathtub. With this patented supporting system you do not only get safely into the bathtub but also safely out of it again.

You can purchase the QP Tub either with or without grip extension. The extension enables broader and thus safer grasping; it efficiently utilizes your energy expenditure and prevents rotating of the body when stepping out.

With the QP Tub you have at your disposition a reliable aid for placing on all smooth and pore-tight surfaces.

Seize the grab bar with both hands before entering or exiting the bathtub and release it only when you sit safely or stand firmly.

Very flexible: The QuattroPower Tub can be attached on the outside or on the inside of your bathtub.




QuattroPower Tub

Optimal. The QuattroPower bathtub systems can be swivelled at three points and thus perfectly adapted to the shape of your bathtub or your bathroom facilities.

Note: The grips are available in two different lengths. When ordering, please take into consideration the dimensions **B** and **C** of your bathtub.

If you have any questions concerning dimensions or possible fields of application of our supporting systems, please ask our team of experts.

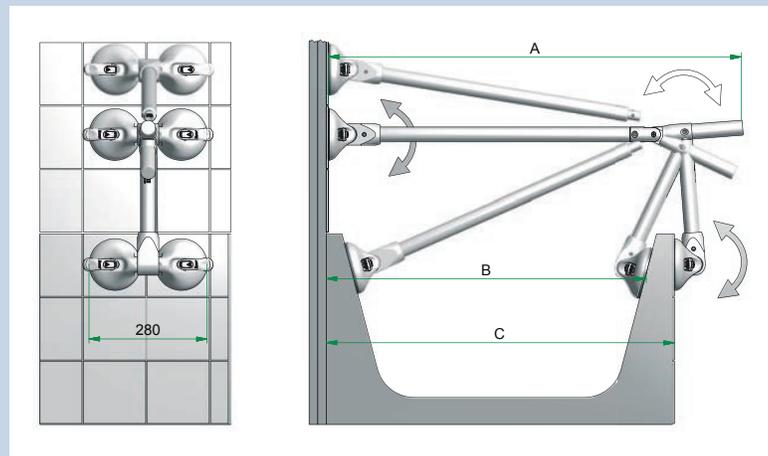

MOBELI	Dimension A	Dimension B	Dimension C	Order No.
QuattroPower Tub (short)	28.8"	24.5"-33.15"	23.4"-28"	R1400263SU
QuattroPower Tub (short) with grip extension	33.15"	24.5"-33.15"	23.4"-28"	R1400264SU
QuattroPower Tub (long)	33.5"	28.8"-38.2"	27.3"-33.5"	R1400261SU
QuattroPower Tub (long) with grip extension	38.2"	28.8"-38.2"	27.3"-33.5"	R1400262SU

CareBag Absorbent Bags

A UNIQUE PATENTED CONCEPT

Hermetic single use bag equipped with a unique super absorbent pad
 Turns urine and soft materials into a gel within seconds
 Designed to safely collect, transport and dispose of body fluids
 Helps prevent the spread of potential pathogens,
 lessening the risk of hospital acquired infections

Also available in



Bedpan Liner

Commode Liner

Vomit Bag

Men's Urinal

Hospitals, Clinics, Long Term Care, Home Health Care, Emergency Preparedness, Hospice, EMTs, Firefighters, School transportation...

aqua Gloves



UNIQUE PRODUCT!

WITHOUT SOAP, RINSING AND DRYING

Pre-moistened with a protective washing lotion dermatologically tested and hypoallergenic

Hospitals,
 Long Term Care Facilities,
 EMTs

1 WASH GLOVE = 1 SHAMPOO!



GelMax Super absorbents

UNIQUE TECHNOLOGY - REDUCES ODORS

Super absorbent powder (available in bulk or water-soluble pouch) and pads which absorb and turn into gel up to 450 times its own weight of fluids and soft materials



Home Care, Hospitals, Labs

Toilevator™ & Toilevator™ Grande

Toilet Base Riser Large Toilet Base Riser





- A clean & safe alternative to standard raised toilet seats
- User's sit on an actual toilet seat
- Eliminates potential balance issues that may occur with other seat risers
- Adds 3.5 inches (9 cm) in height for comfort and ease
- Discreet & durable; fits under existing toilet base
- The Kits come complete with all parts included
- DIY installation is easy.
- No special cleaning or maintenance required
- Fits standard round or elongated toilet bases
- Patented and award winning design

For information on any of these products from HartMobility please contact Clarke Health Care Products.



NEW The all new **Toilvator™ GRANDE** fits most of the new, larger base toilets such as American Standard Boulevard, & Antiquity series, Kohler Cimarron, Bancroft and Rialto series, Crane Titan series, Eljer Titan Pro series, and many more.

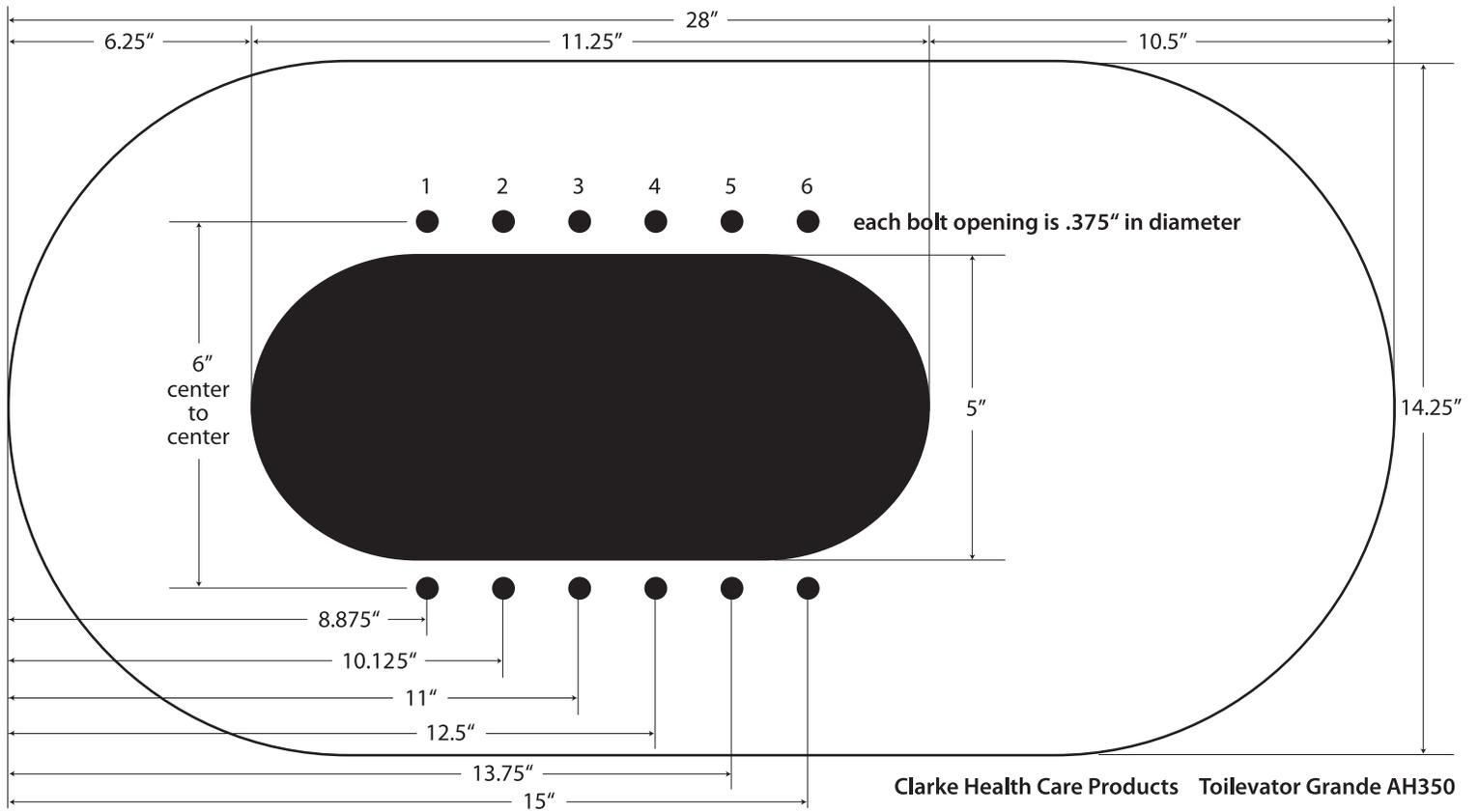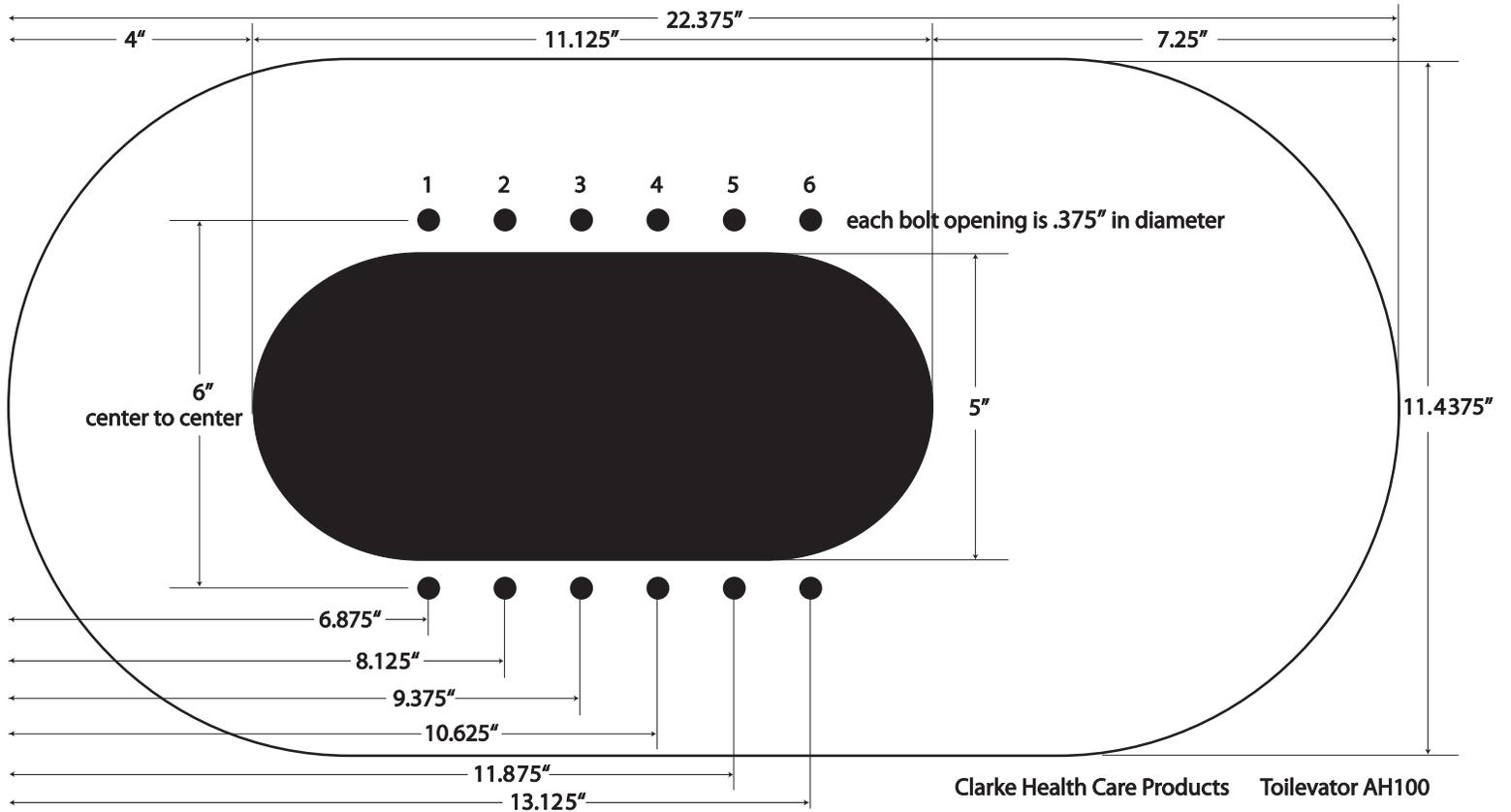

www.clarkehealthcare.com
888-347-4537



The LEGACY

Dolomite developed the Legacy to accommodate a wide variety of users with the most standard features for the best value. The design incorporates these features: flip-up seat, folding frame, adjustable handles, and wide wheels. With choice of color and an extensive accessory line, the Legacy looks smart strolling down the sidewalk.

Frame

The strongest frame on the market, Dolomite welds the steel frame for greater durability, no screws to loosen at joints. Legacy is available in two colors: Metallic Blue or Candy Apple Red. There are three frame sizes to accommodate individuals of all sizes. Dolomite frames are easily folded for storage or transportation. The handles are height-adjustable in 4 increments.

Braking

Legacy comes with the patented Dolomite Automatic Braking System with adjustable tension for varied grip strength. Brake handles have wide flange and large loops for comfortable hand positioning. When parked and using the seat, push down the handle to set the parking brakes. The brake mechanism is located within fork for protection from damage. Accessory slow down brake provides adjustable tension for hilly terrain or to even out direction of travel for stroke patients.


Standard Features

Legacy supports 350 pounds, which is 75 pounds more than the majority of the other rolling walkers. Other standard features are:

- Removable basket with carry handle
- Flip-up seat for greater reach or gait training
- Seat handle releases to fold walker
- Foam grips for comfortable and secure grip
- Flexible and adjustable back strap
- 8" indoor/outdoor tires

Warranty

Dolomite provides a two year warranty against manufacturer's defects. This excludes grips, cables and tires.

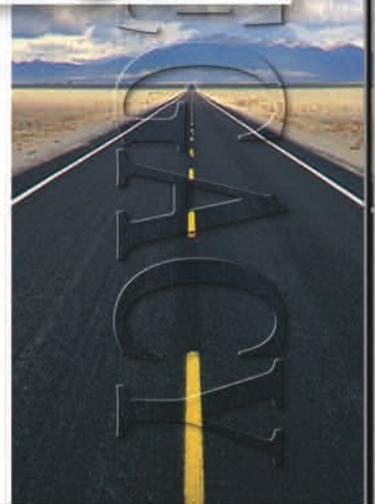

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com



The LEGACY

Accessories



OXYGEN TANK HOLDER



ONE-HAND BRAKE



SLOW DOWN BRAKE



CROSS BAR WITH BRAKE



FOOD & BEVERAGE TRAY



CANE HOLDER



TALL HANDLES



UNDER SEAT BAG



BACKSTRAP PAD,
SEAT CUSHION


Accessories can enhance the use of your walker – the Tray gives you a platform to carry food, the Oxygen Tank Holder will safely support a B, C or D tank. A Cane Holder can secure a cane or crutch to the side of the walker, Slow Down Brakes have a tension control to allow individual settings for each wheel. One-Hand Brakes can control your braking from one handle. Under Seat Bag gives you a place to stash items out of sight. Tall handles add 3” to the handle height. Cross bar combines more hand hold area with one hand braking.

SPECIFICATIONS	SUPER LOW	LOW	STANDARD
Dist. Between Rear Wheels	19"	19"	19"
Minimum & Maximum Handle Height	28 - 31"	31 - 34"	33 - 37"
Minimum & Maximum Tall Handle Heights	31 - 35"	33 - 37"	35 - 39"
Seat Height	19"	22"	24"
Overall Width/Length	23"/23"	23"/26"	23"/26"
Width Between Handles	17"	17"	17"
Folded Dimensions	36" Length x 23" Width x 12" Height		
Walker Weight (lbs.)	16	17	17
Weight Capacity (lbs.)	350	350	350
Product Number: Blue	D12054B	D12052B	D12050B
Product Number: Red	D12054R	D12052R	D12050R

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

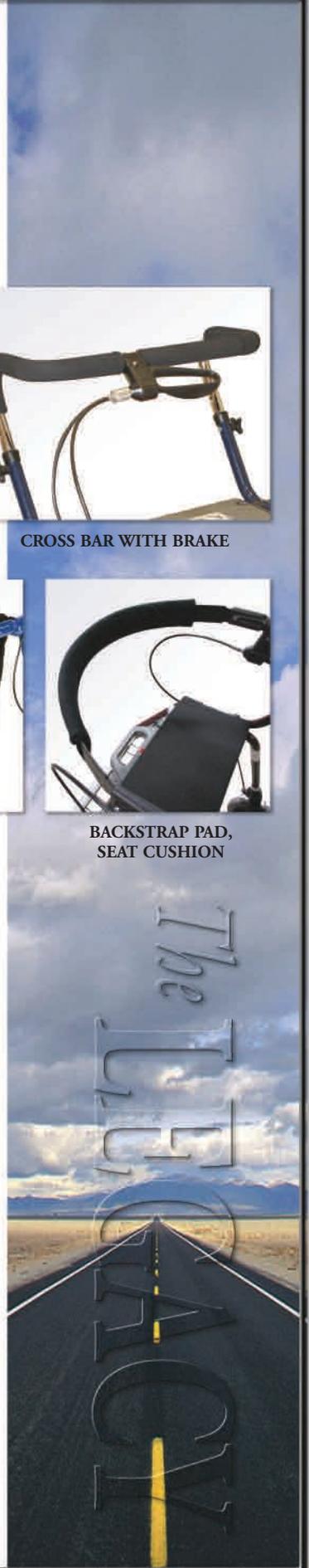



The MAXI+

Maxi+ is designed for the heavier, larger, or taller user. The wider frame provides more seat room and a larger walking area. Dolomite has the strongest frames on the market; Maxi+ is rated for a 440 lb. capacity. Choose Maxi+ for a stable walking experience and travel safely in style.

Components

Maxi+ is made from high tensile steel with welded joints for strength and low maintenance. Maxi+ comes in Metallic Blue epoxy-based powder coat finish. Dolomite Maxi+ comes in two frame heights; see chart for specs. The frames fold front to back for storage. A safety latch prevents accidental folding when using the seat handle.

Maxi+ includes a wide flip-up seat for more walking room. Anatomical grips have a larger surface, which reduces hand fatigue for arthritic users. Flexible back strap is adjustable. Large 8" wheels roll easily indoors or outdoors. Red Curb climber on the rear wheel assists when lifting the walker over a curb or step.

Braking

Maxi+ comes with the patented DABS braking system, with adjustable tension for varied grip strength. Brake handles have a wide flange and large loops for comfortable hand positioning. The brake mechanism is located within the fork for protection. When parked and using the seat, push down the handles to set the parking brakes.

Features

- Extra wide, 20.5" between handles
- Large flip-up seat
- 440 pound capacity
- Perfect for the taller person, 25.5" seat height
- Large 8" wheels
- Handles adjust 7" in 1" increments
- Large basket for personal items
- Adjustable back support

Warranty


Dolomite provides a two year warranty against manufacturer's defects. This excludes grips, cables and tires.



CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

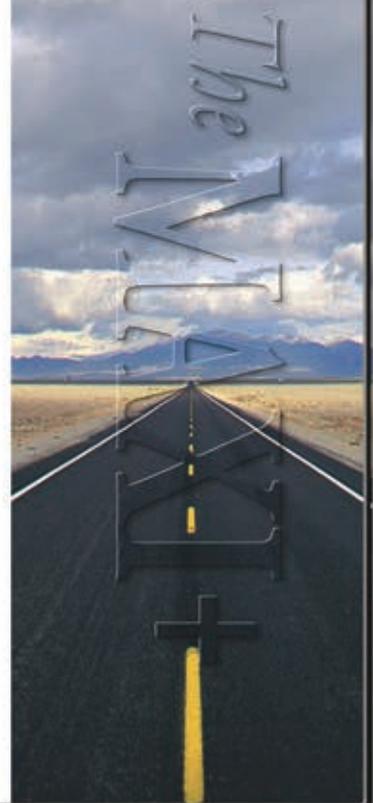



The MAXI+

Accessories



TRAY



OXYGEN TANK HOLDER



CANE HOLDER



SLOW DOWN BRAKE



ONE HAND BRAKE



UNDER SEAT BAG

Accessories can enhance the use of your walker – the Tray gives you a platform to carry food, the Oxygen Tank Holder will safely support a B, C or D tank. A Cane Holder can secure a cane or crutch to the side of the walker, Slow Down Brakes have a tension control to allow individual settings for each wheel. One-Hand Brakes can control your braking from one handle. Under Seat Bag gives you a place to stash items out of sight.


Specifications:

SPECIFICATIONS	LOW	STANDARD
Dist. Between Rear Wheels	22"	22"
Minimum Handle Height	28.5"	32.5"
Maximum Handle Height	35"	39.5"
Seat Height from Floor	22"	25.5"
Seat Dimensions	8.5" Width x 17" Height	
Overall Width	27.5"	27.5"
Overall Length	26"	29"
Width Between Handles	20.5"	20.5"
Folded Dimensions	33" Length x 27.5" Width x 11" Height	
Walker Weight (lbs.)	20.5	21
Weight Capacity (lbs.)	440	440
Product Number	D12122	D12120

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

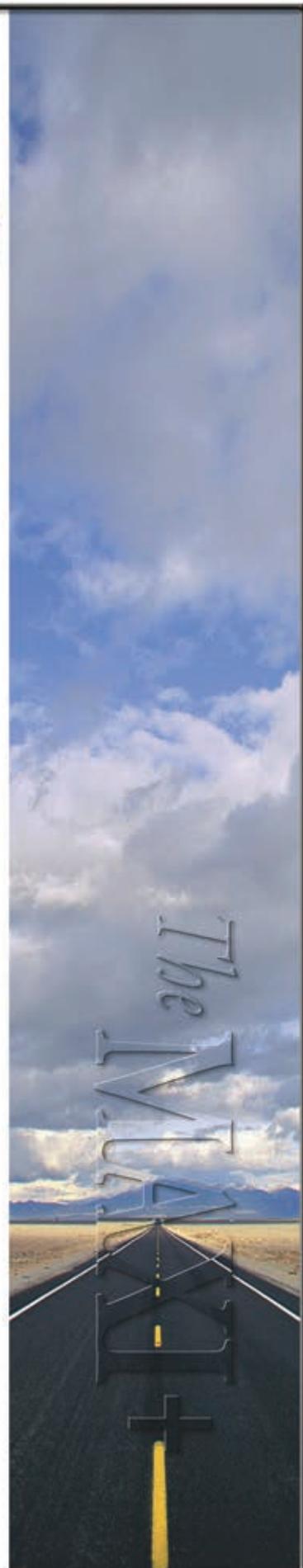



The SYMPHONY

Symphony is a elegant orchestration of quality and durability and is the finest rollator on the market today. The handles have stepless adjustment and anatomic hand grips. These grips require less grip strength and reduce hand fatigue especially for arthritic hands. **Symphony** has a flip-up seat for greater reach or for gait training.

Frame

The lightweight aluminum frame is welded for greater durability. No screws to loosen at the joints. **Symphony** is the strongest aluminum frame available today. Grey anodized finish. **Symphony** is available in 3 frame heights. Easy to fold frame has a safety latch to prevent accidental folding. Aluminum frames are suitable for corrosive environments, such as coastal areas.

Braking

Symphony comes with the patented **Dolomite** Automatic Braking System with adjustable tension for varied grip strength. Brake handles have wide flange and large loops for comfortable hand positioning. When parked and using the seat, push down the handles to set the parking brakes. The brake mechanism is located within the fork for protection. Accessory slow down brake provides adjustable tension for hilly terrain or to even out the direction of travel for stroke patients.

Standard Features

- Flip-up seat for gait training and greater reach
- Curb climber on rear wheel assists when lifting front of walker
- Anatomic grips reduce hand fatigue; can be rotated inward for a more natural position
- Adjustable back support
- Basket for shopping
- 8" indoor/outdoor wheels
- Stepless handle height adjustment, adjustment marks printed on handle

Warranty

Dolomite provides a two year warranty against manufacturer's defects. This excludes grips, cables and tires.



CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

The SYMPHONY



The SYMPHONY

Accessories



OXYGEN TANK HOLDER



ONE-HAND BRAKE



SLOW DOWN BRAKE



UNDER SEAT BAG



FOOD & BEVERAGE TRAY



CANE HOLDER



CROSSBAR



HEMIPLEGIC HANDLE


Accessories can enhance the use of your walker – the Tray gives you a platform to carry food, the Oxygen Tank Holder will safely support a B, C or D tank. A Cane Holder can secure a cane or crutch to the side of the walker, Slow Down Brakes have a tension control to allow individual settings for each wheel. One-Hand Brakes can control your braking from one handle. The Hemiplegic Handle allows you to use the body to push for momentum while braking is done with one hand. Under Seat Bag gives you a place to stash items out of sight. Cross bar provides one hand braking with a large push bar.

FOLDED SPECIFICATIONS	SUPER LOW	LOW	STANDARD
Dist. Between Rear Wheels	18.5"	18.5"	18.5"
Minimum Handle Height	25"	27.5"	30"
Maximum Handle Height	30"	33"	36"
Seat Height	18"	21"	24"
Overall Width	23"	23"	23"
Overall Length	22"	24"	27"
Width Between Handles	17"	17"	17"
Walker Weight (lbs.)	15	15	15
Capacity (lbs.)	275	275	275
Product Number	D12162	D12161	D12160

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

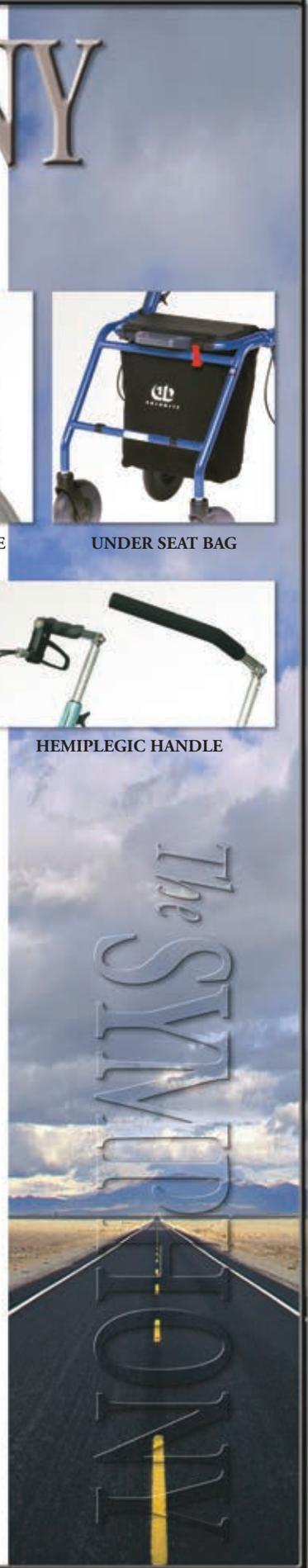



The JAZZ

Linking Style with Function

The **Jazz** is an innovative, lightweight rollator designed to offer excellent support and stability. The new stylish rollator includes numerous smart features making it the ideal walking companion. It's ergonomic design allows a natural walking position, within the rollator.



Features

- Padded seat for comfort
- Handle height memory, remembers the preset height.
- Design provides more walking space within the rollator
- Cloth basket and back strap included

Components

The newest **Jazz** is constructed with a welded aluminum frame for weight savings. Welded joints do not need to be maintained like a frame assembled with hardware. A safety latch locks the **Jazz** in the open, walking position. Lift up on the seat strap to easily fold the **Jazz**. The side to side fold allows it to stand by itself. The small profile can be easily stored in your home. Once set, the handle memory feature allows you to lower the handles for storage, then raise them back to the preset position without any additional adjustment. The handles have a pressure-relieving, anatomic grip.

Squeeze the **Jazz** brake handles to slow down and push down to set the parking brake for sitting. Brake cables are fully integrated into the frame for a clean look and to prevent catching on objects while walking. Eight inch wheels are suitable for indoor and outdoor travel. A molded curb climber is mounted on the right rear wheel. Placing your foot on this anchors the wheel to lift the front of your walker over a curb or obstacle. Can be moved to left wheel.

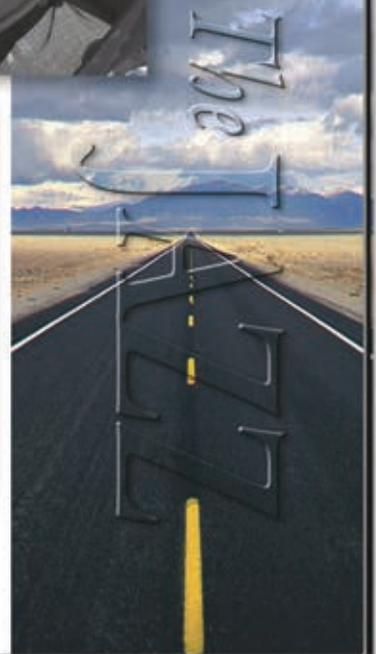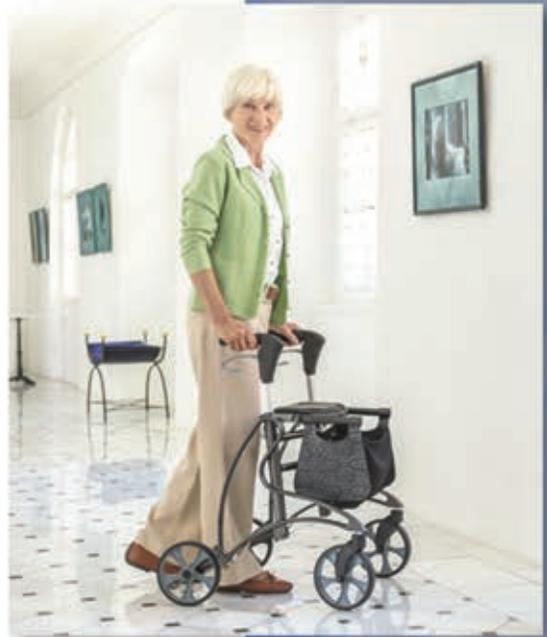

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com



The JAZZ Accessories



BIG BASKET, BLACK



BIG BASKET, RED



BACK STRAP PAD

Accessories

- Tray
- Big Basket, Black
- Big Basket, Red
- Back Strap pad
- Cane Holder
- Slow Down Brakes



TRAY



CANE HOLDER

Warranty

Dolomite provides a 2 year warranty against manufacturer's defects. This excludes grips, cables and tires.


Specifications

SPECIFICATIONS	STANDARD 610	LOW 510
Basket capacity	10 lbs	10 lbs
Wheel	8" D x 1 1/4" W	8" D x 1 1/4" W
Capacity	330 lbs	330 lbs
Weight of walker	17 lbs	16.5 lbs
Overall width	23.2"	23.2"
Overall length	26.9"	26.9"
Width between wheels	19.9"	19.9"
Overall height	29.25" – 39"	26.35" – 31.2"
Seat height	23.8"	19.9"
Seat	14 1/2" L x 8" W	14 1/2" L x 8" W
Width between handles	17.5"	17.5"
Folded width	8.3"	8.3"
Product Number	D1559202	D1559201

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

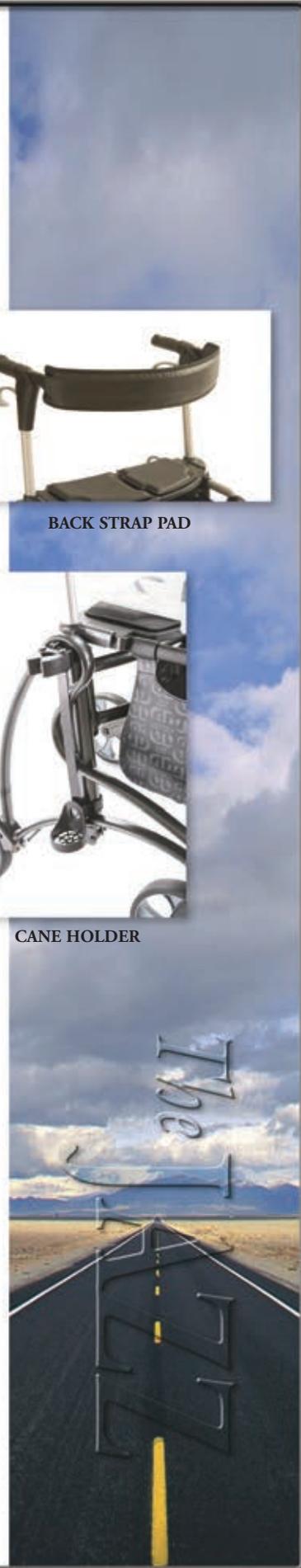



The FUTURA Tall

FuturaTall has been designed to accommodate persons taller than 6'2" who need walking assistance. It has the tallest handles in the market on a standard rollator, at 43.5", and a seat height of 27.5". Incorporating all the features you have come to expect in a Dolomite rollator, **FuturaTall** fills the specific needs of a taller population.

Frame

The strongest frames on the market, **Dolomite** welds the steel frames for greater durability – no weak mechanical joints. **FuturaTall** comes in a sophisticated Black-Grey color. The frame folds for storage and transportation. Handles are infinitely adjustable within a six-inch range. **FuturaTall** has a 330 lb. capacity. Flip-up seat gives greater walking clearance.

Braking

Dolomite patented braking system has adjustable tension for varied hand strength. Brake handles have a wide flange and large loops for hand positioning. A stainless steel brake located within the fork gives you years of service. Push downward to set a parking brake while seated. Accessory slow down brakes provide additional tension on the wheels for hilly terrain or for assistance while gait training.

Standard Features

- Flip-up seat for greater reach
- Seat handle release to fold
- Magnetic seat latch secures against accidental closure
- Foam covered handles for comfortable grip
- Flexible and adjustable back strap
- Removable basket with carry handle
- 8" indoor-outdoor wheels

Warranty

Dolomite provides a two year warranty against manufacturer's defects. This excludes grips, cables and tires.



Fully adjustable handles

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com



The FUTURA Tall Accessories



SLOW DOWN BRAKE



UNDER SEAT BAG



OXYGEN TANK HOLDER



TRAY

Accessories

The extensive line of accessories enhances the use and comfort of your walker. For instance, a molded plastic tray fits onto the seat to carry food or beverage. The oxygen tank holder will replace the basket and securely carries a B, C or D size tank. You can secure a crutch or cane on the side of your walker with the 2-piece cane holder. Slow down brakes have a separate tension control to create drag on your wheels, slowing them for gait training, stroke patients or useful on hilly terrain. When not using a basket, the under seat bag gives you a place to hold personal items out of site. This bag can be used in conjunction with the oxygen tank holder. Anatomic grips are a formed rubber grip shaped to the palm of your hand. The one-hand brake makes it possible to activate both brakes from one side.



CANE HOLDER



ONE HAND BRAKE



CURB CLIMBER

Specifications

DESCRIPTION	SPECIFICATION
Distance between rear wheels	21.5 in.
Minimum handle height	37.75 in.
Maximum handle height	43.5 in.
Seat height	27.5 in.
Seat L x W	14 x 7 in.
Width between handles	16.75 in.
Overall width	25.5 in.
Overall length	30 in.
Folded dimensions	44.25 in. (L) x 25.5 in. (W) x 12.5 in. (H)
Weight without basket	18 lbs.
Capacity	330 lbs.
Product Number	D12076T

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com





The SOPRANO

Amid a chorus of *oohs* and *ahhs*, the **Soprano** has joined the line of beautifully engineered **Dolomite** walkers. **Soprano**, with its larger wheelbase, provides a smooth ride over life's bumps. Completely redesigned with new frame, seat, handles and cable-less braking system.

Frame

Soprano is based on a welded steel frame. The frame is designed to eliminate the rear forks, giving you more walking space between back wheels. A larger seat flips up to give you more walking room. The front wire platform can carry up to 26 pounds of baggage.

Braking


The brake handles and wedge shaped hand pad offer a secure and comfortable grip. Cable-less braking system uses stainless steel rods encased in the frame to actuate the brake. The brake shoe is designed to catch dirt and stones before they get stuck in the brake pads, adding longevity to the wheels and brakes. Infinitely height adjustable handles also rotate for better wrist positioning.

Standard Features


- Welded steel frame
- 10" wheels
- Wire shelf carries your bulky items
- Handles rotate and are height adjustable
- Champagne frame color
- Basket with carry handle

Warranty

Dolomite provides a two year warranty against manufacturer's defects. This excludes grips, cables and tires.

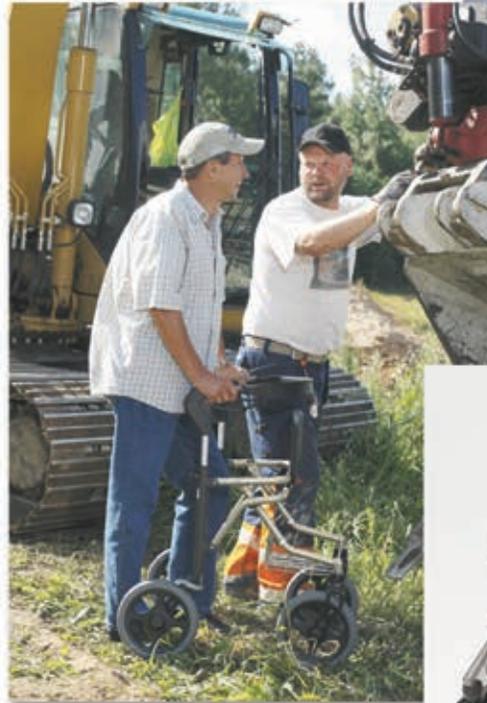

Folded position

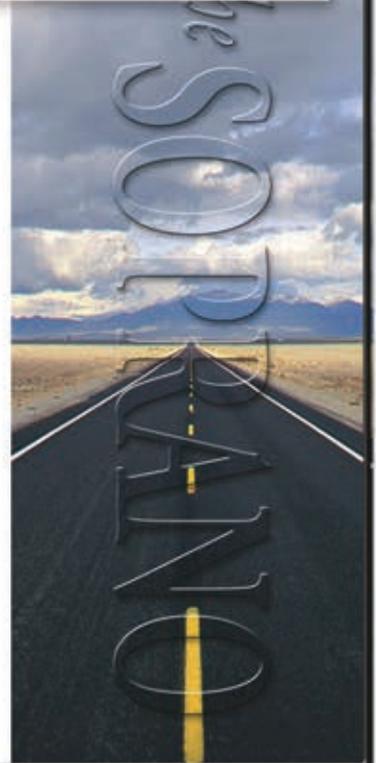

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com



The SOPRANO

Accessories



CANE HOLDER



WHEEL COVER



CURB CLIMBER



TRAY



O2 TANK HOLDER



SLOW DOWN BRAKES

The following accessories are offered for your **Soprano**. Cane holder attaches to the frame to hold a cane or crutch. The wheel cover keeps dirt or mud from splattering on your clothing. Use the curb climber to assist when tilting the walker to climb over curbs. Food and beverage tray sits on the seat, will also hold small items. Safely carry a B, C or D size tank with the Oxygen Tank Holder. Slow-down brakes add tension to the wheels so the walker won't roll away, useful for gait training or hilly terrain. The tension is self-adjustable and can be secured in position.

SPECIFICATIONS	SOPRANO
Inside width between rear wheels	20"
Minimum handle height	30.5"
Maximum handle height	37"
Seat height	23.5"
Width between handles at seat	17"
Overall length of walker, front to back	33.5"
Overall width	25.5"
Folded dimensions	11.5"H x 25.5" W x 34.5" L
Walker capacity	355 lbs.
Walker weight	21.4 lbs.
Product number	D12030

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com





The ALPHA

Dolomite's Alpha walkers offer walking support for people with arthritis, neurological disorders, Parkinsons, stroke survivors or people with limitations in wrist or hand movement. A variety of adjustments can be made for arm supported walking. You can use forearms to propel the walker when there is limited hand strength. Arm height adjusts over 9" and the rotation of the arms adjusts 20°.

Components

Dolomite uses a welded steel frame, no hardware to loosen at joints. Alpha frame has been redesigned wider and deeper to accommodate the higher center of gravity caused by the arm supports. Alpha has 8" high friction, indoor/outdoor tires. Rear of padded forearm supports is recessed to relieve pressure on the elbow. An Allen key is included for adjustments to the arm supports. Step-less arm height adjustment has measurement marks printed on tube to duplicate settings. Flip-up seat gives the user the ability to "walk inside" the frame for better gait training.

Braking

Alpha Advanced brakes use a joystick style grip. Brake handles can be adjusted so that little force or leverage is needed. A ball joint on allows you to change the direction of brake operation. Alpha Basic has a loop brake handle.

Standard Features

- Folding frame
- Flip-up seat for gait training and greater reach
- Curb climber on rear wheel assists when lifting front of walker
- Removable basket with carry handle
- Magnetic security latch prevents accidental closure
- Attached Allen key for arm adjustments

Warranty

Dolomite provides a two year warranty against manufacturer's defects. This excludes grips, cables and tires.



BASIC



ADVANCED



CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com



The ALPHA

Accessories



OXYGEN TANK HOLDER



ONE-HAND BRAKE



SLOW DOWN BRAKE



FOOD & BEVERAGE TRAY



CANE HOLDER



PLATFORM SUPPORT

Dolomite has extensive accessory offerings. Accessories can enhance the use of your walker; for instance, the Tray gives you a platform to carry food and the Oxygen Tank Holder will safely support a C or D tank. You can secure a cane or crutch to the side of the walker with the Cane Holder. Slow Down Brakes have a tension control to allow individual settings for each wheel. One-Hand Brakes control standard brakes from one handle. Platform Support replaces the arm troughs; large padded surface to lean onto. Use with either brake system.

SPECIFICATIONS	ALPHA ADVANCED / BASIC
Overall length	27.5"
Overall width	25.5"
Distance between handles	12" - 14"
Distance between rear wheels	20"
Seat height	24.5"
Forearm pad to floor	33.9" - 46.4"
Capacity	275 lbs
Walker weight	28 lbs.
Product Number: Advanced	D12220
Product Number: Basic	D12230

CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com



DOLOMITE | *STEP UP*



The rehab optimizer

DOLOMITE STEP UP is a new modern platform support. Made to give rehabilitation to several patients every day. Adjust to optimal ergonomic position by a few simple handgrips – quick and easy!

DOLOMITE STEP UP can be manually adjusted or hydraulically. Assembly of accessories is tool free. Are you ready to step up?





Made for training

The base frame is designed to enhance premium gait area. The raised frame in front allows longer step. The width in side the frame generous.

Ergonomic adjustment of arm support

Arm pads are in soft material for pressure relief. Securely assembled with microfilm. Quick and simple to remove for cleaning or replacement.

Transfer and access

Well thought out for transfer situations. From bed, the lowered rear-end helps reversing under the bed. When

transferring from wheel chairs benefit of the generous frame size and, the width adjustment.

Features and Options



By releasing the central knob, the both arm pads are simultaneously adjusted in width. The position of the arm supports remains in an angle of 25 degree from the shoulder points - the ideal position for keeping upper strength.

Accessories and options

Quick and easy assembly in track or prepared recess.

**Practical basket**

Keeping belongings within easy to reach range.
Art.no **D1531390**

**Tray**

for small things.
Art.no **D1531391**

**Elbow stop**

Securing elbows from falling out of position.
Art.no **D1531448**

**Swivel caster with directional brakes**

Helps focusing on walking instead of direction.
Art.no **D1531408**

**Brake system**

For the user capable of operating brake function.

Hand-brake kit

Art.no **D1537656**

One-hand brake

Art. no Right
Art. no Left

D1539513
D1531450

**Oxygen tank holder**

Secure mounting of oxygen tube on the walker.

Art.no **D1531406**



The sling gives security and make it possible for the user to start the walk training earlier. No tools needed. Easy to fit and adjust, quick release.
Art. no **D1536846**

**Stand assist incl. slings**

For verticalization.
Art.no **D1531449** (standing plate + knee support + sling)

Standing plate

Art.no. **D1535330**

**IV holder**

Separate IV-stand makes life complicated. Not anymore.
Art.no **D1531407**

**Ball grips**

when having difficulties to close the hand.
Art.no **D1539515**

Small caster 100 mm diameter.
Art. no **D1536047**

Others

Models



Manual

Pull out the pin move to required position release the pin.

Pull up Hydraulic

Pull up lever to adjust to required height.

Technical data



22.6"/26"/29.6"



27.3"/30.8"/34.3"



34" - 51.5"



13.2" - 18.3"



24" - 29.25"

Step Up Manual

Step Up Hydraulic



31.5"



35.2 lb



330 lb



5"



8.3"

Step Up Manual

Step Up Hydraulic

31.5"

39.6 lb

330 lb

5"

8.3"

Colour



RAL 1035

Article code

Step Up Manual D1521677
Step Up Hydraulic D1521678



www.clarkehealthcare.com
 888-347-4537





INVACARE® WALKERS

GSA Schedule



With the **Invacare® Walker**, providers, clinicians and consumers gain the stability they need at an economical price. Invacare Walkers offer a wide, deep frame with a large number of height adjustments. The walker has a lower side brace for added stability and is stable, lightweight and easy to lift and maneuver. Walkers include PVC handgrips that provide comfort and long-lasting wear and anti-rattle silencers for quiet operation.

Invacare® IClass™ Paddle Walker Features

Lightweight, wide, deep aluminum frame

Composite lower side brace stiffens walker frame and adds stability

Dual Blue-Release mechanisms provide both visual and audible "locked" cues

Folding mechanisms are easy to reach from a seated position

Lighter and easier to lift and maneuver

PVS handgrips wear longer and provide greater comfort

Large number of height adjustments

Anti-rattle "silencers" for quiet operation

Walkers fold easily for transport or storage

HCPCS Code:

E0135 Adult/Jr./Single

E0143 Wheeled Adult/Jr.



1. Dual-Release Paddle button



6291-A: Dual-Release Paddle Adult Walker



6291-5F: 629171, Dual-Release Adult Paddle Walker with 5" Wheels

Walker Specifications



Model #	Width Inside Hand Grips	Width Inside Base Legs	Depth Opened At Base	Depth Folded	Height Adjustment	Product Weight	Weight Capacity
6291-1	17"	21.5"	17.5"	4"	30.375"-37.375"	5 lb. 5 oz.	300 lb.
6291-5F**	17"	21.5"	17.5"	4.5"	33"-39"	5 lb. 6 oz.	300 lb.
629171	17"	21.5"	17.5"	4.5"	33"-39"	5 lb. 6 oz.	300 lb.




	Patient Height Range*	Units Per Carton	Limited Warranty
6291-1	5'3"-6'4"	1	5 Years
6291-5F**	5'6"-6'6"	4	5 Years
629171	5'6"-6'6"	1	5 Years

* All patient height ranges are approximate. Each individual should always consult with his/her physician or therapist to determine proper adjustment and usage.

** With 13" rear leg extensions.

"F" designation in Model no. 5" fixed wheels with glide tips.

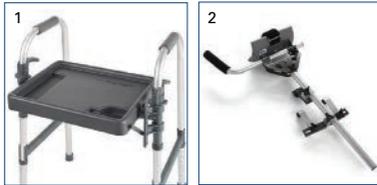

1. **Walker Tray**
Model # 6007
Compatible with all walkers
2. **Platform Attachment**
Model # 6027
Compatible with all walkers



1. **Glide Tips**
Model # 6280
Packaged 10 per carton
Durable construction for long wear



1. **5" Single Fixed Wheel Attachments with Rear Glide Tips**
Model # 6271
Use 13" rear leg extensions for maximum adjustability



1. **Tall Leg Extensions**
Model # 6275
Adjusts to accommodate persons up to 6' 6" tall
2. **13" Rear Leg Extensions with Glide Tips**
Model # 6266
Use with swivel or fixed wheels

Accessory Specifications



Model #

Patient Height
Range Adult
min. – max.

Patient Height
Range Junior
min. – max.

Limited
Warranty

Fixed-wheel

6271 5" single fixed-wheels w/rear glide tips	33" – 39"	28" – 33"	5 Years
--	-----------	-----------	---------

Leg extensions

6266 13" rear leg extensions w/glide tips	33" – 40"	29" – 35"	5 Years
--	-----------	-----------	---------

6275 tall front and rear leg extensions	33" – 40"	28" – 35"	5 Years
--	-----------	-----------	---------

Clarke Health Care
www.clarkehealthcare.com



©2012
Invacare Corporation. All rights reserved. Trademarks are identified by the symbols TM and [®]. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Lotus is a registered trademark of Inovo, Inc. Specifications are subject to change without notification.

Alpha Advanced and Alpha Basic



The Alpha models have forearm supports with brake handles at the end of the pads. The Advanced is a joy-stick style handle which is push and pulled. It can be rotated on the ball joint mount to also be pushed side to side for braking and setting the park.

The Basic handle is the Dolomite looped brake handle mounted vertically. You squeeze to brake and push out to set the park.



Alpha models can also be built with one hand brakes on either side, using either style brake handle.



A large pad can be substituted for the individual troughs. Similar to a 3 sided tray, it permits the user to lean in at any point and reach across to the brake handles in front.



Walker base models

We use the Legacy, Maxi or Symphony as our bases to build these custom rollators. The Legacy and Symphony come in 3 seat heights. These models allow us to configure a walker with forearm pads on a smaller scale. The seat height of the super low models is only 18" and 19". The Legacy also has the option for us to use the tall handles to gain 3" additional handle height.

The Maxi+ has a capacity of 440 lbs., a wider seat and comes in 2 seat heights.



To order a custom Dolomite, provide some details about your client, weight, height, ability or special need. We'll design a system and send you a quote.

We would prefer to install these accessories and ship you a completed walker. But with some instruction these can be added to an existing walker.

Contact us:

Phone: 888-347-4537

fax 724-695-2922

email info@clarkehealthcare.com,



Clarke Health Care Products

Let us customize a Dolomite walker for your client's special needs.

Using the Dolomite Legacy, Symphony or Maxi as our base, we can add accessories or modify handles and brakes to accommodate their ability to use a rollator.



Clarke Health Care Products

7830 Steubenville Pike

Oakdale PA 15071

Phone: 888-347-4537

Email: info@clarkehealthcare.com

www.clarkehealthcare.com

One Hand Brakes



For amputees, stroke survivors, one hand brakes allows control of the braking action from one hand. The two brake

cables will exit the brake handle, with one cable crossing over to the other brake pad. The other handle can be a forearm trough or only a handle to grasp.

Crossbar



Another one hand brake option is the Crossbar, this is a tubular bar that replaces the handles and back strap. It resembles the push bar of a shopping cart with a brake handle in the center. The padded bar provides larger area to grasp and push the rollator. The looped brake handle in the center can be a right or left one-hand brake system.

Forearm Troughs



Forearm pads, or troughs allow the client to lean on their forearms with brake control at the end of the pads. The brakes can be on mounted on both sides or controlled from one side. The pads are lightly cushioned and have some adjustment in angle and position in front of the body.

Note that using the forearm pads, you restrict use of the seat. Depending on position and height, the pads may be in the way to comfortably sit on the seat.

HP Handles



HP handles or hemiplegic handles are mounted on the Symphony model. The one hand braking system uses a padded bar on the non-braking side. This bar can be adjusted in height and angle to the body, so that you can push with your body while steering and braking from the other side.

Slow Down Brakes



If the walker seems to be moving too easily or getting away from the client, add Slow Down brakes. They do not affect the braking action, but only add tension or drag on the wheels. The thumb screw will tighten a nylon wheel onto the back tire. The left and right sides can be individually adjusted. For stroke survivors who may be stronger on one side, tighten that side so they walk in a straight line. For gait training, tighten both to slow the tire rotation. You will still have all the braking action from the brake handles. There is a hidden nut under the thumbscrew to lock in your adjustment, if needed.





Custom Dolomite Rollator Worksheet for quote, return to info@clarkehealthcare.com or fax: 724-695-2922

1. Walker	Choose one model, size or color	4'10"-5'2" suggested height	5'2"-5'5" suggested height	5'4"-6'2" suggested height
	Legacy , <330 lbs cap. <div style="text-align: right;">Blue</div>	<input type="checkbox"/> Legacy 450 SuperLow Seat height 19", handle 28-31"	<input type="checkbox"/> Legacy 520 Low Seat height 22", handle 31-34"	<input type="checkbox"/> Legacy 600 Standard Seat height 24", handle 33-37"
	Symphony, <275 lbs cap. Silver frame	<input type="checkbox"/> Sym. 450 SuperLow Seat height 19", handle 25-30"	<input type="checkbox"/> Sym. 520 Low Seat height 22", handle 27.5-33"	<input type="checkbox"/> Sym. 600 Standard Seat height 24", handle 30-36"
	Maxi+, <452 lbs cap., Blue frame		<input type="checkbox"/> Maxi+ 550 Low Seat height 22", handle 28.5-35"	<input type="checkbox"/> Maxi+ 650 Standard Seat height 25.5", handle 32-38"

2. Arm troughs or platform , does not include brake handles, choose brake below	<input type="checkbox"/> One pad left	<input type="checkbox"/> One pad right	<input type="checkbox"/> 2 Pads, left & right side
Option 1 Forearm pad, adds 3" height	<i>If using pad without a brake, we will add a grip at end of pad</i>		
Option 2 Platform, adds 4" height	<input type="checkbox"/> Full platform, brakes mount at front edge, 27"w x 15.5"h		

3. Brake handle position and style	One hand, operates brakes on both wheels	Braking both hands
A. Choose brake position	<input type="checkbox"/> Left mount one hand brake	<input type="checkbox"/> Right mount, one hand brake
B. Choose style of brake handle and action	<input type="checkbox"/> Joystick style, push/pull	<input type="checkbox"/> Looped handle, squeeze/push
<input type="checkbox"/> Brakes on both sides		

Alternate one hand braking systems, will not work with arm pads above			
Crossbar, front push bar w/ one hand brake mounted center, replaces back strap, <i>not available for Maxi+</i>	<input type="checkbox"/> Left brake mount	<input type="checkbox"/> Right brake mount	
HP Handle, one hand brake w/ angle adjustable bar on other side, <i>only available on Symphony, Left or right in prod # is based on HP side.</i>	<input type="checkbox"/> Left brake, right HP	<input type="checkbox"/> Right brake, left HP	


Accessories	D12533 Tray	D12682 O2 Tank Holder	
	D13041 Slow Down Brakes	P12519 Seat Pad,	
	D12606 Cane holder	D12523 Under seat bag	
	D13235 Tall Handles for Legacy, adds 3" height		
	D13235M Tall Handles Maxi, adds 4" height		

Ordering company _____
Return quote to _____
Email _____@_____
Ph: _____ Zipcode for delivery _____

Mark For _____ User weight _____ User height _____ Floor to Elbow _____

This walker will be built in our shop. It cannot be shipped totally assembled. The handles will need to be inserted and adjusted for height and length. It is not recommended for drop ship delivery unless a technician can be present to make the adjustments to fit the client. Arm troughs may block use of seat for sitting, platform will block use of seat. Note, the part numbers in the chart above do not indicate all components necessary to complete the walker assembly. Quote will show all necessary pieces. Non-returnable.

Custom options

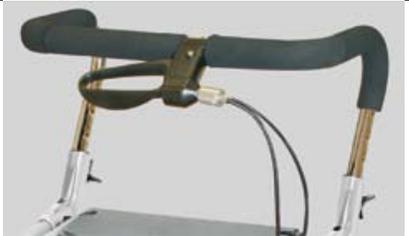Forearm trough or Pad options		
		
Forearm troughs/pads. Adds 3" to handle height. Choose brake handle style	Full platform showing joystick brake handles. Adds 3" to handle height. Choose brake handle style	
Brake handle options used with arm pads or platform		
		
Looped brake handle, used with arm pad or platform, showing one hand braking, can also be mounted for 2 hand braking. Squeeze and push brake action	Joystick style braking used with arm pad or platform, can be mounted as one hand or two hand braking, Push/pull brake action.	
One hand braking		
		
Cross bar, replaces back strap, brake handle mounted center. Left or right one hand braking	One hand brake, 2 cables from one side operate brakes at both wheels. Can be used on Legacy, Maxi or Symphony without pads or platform.	HP handle, one hand brake with padded bar on other side. Replaces backstrap. Push with trunk or forearm. Only used with Symphony

The DecPac Portable Ramp

DecPac Fiberglass Ramps

Wherever your entry problem, front porch, door threshold, curb or minivan, DecPac folding ramps are a light-weight solution. No assembly needed, ready to use right out of the box, simply unfold. Great for wheelchairs, walkers, strollers or anyone with limited mobility. No permanent assembly or installation, quick to setup or remove. Perfect for renters or travelers.

Features

DecPac ramps are lightweight and quiet, will silently ride in the back of your vehicle. All sizes of ramps have a built-in handle for carrying. The ends of the ramp have a ribbed grip-edge so they stay in place. Accordion-style fold makes the ramps easy to fold and compact so that they slide under a minivan seat or stand in a closet.




Components

Constructed of fiberglass-laminate with a rigid foam core, DecPac ramps are weather resistant. The deck has non-slip surface and bright safety tape to mark the edge. Personal size ramps have 3 panels. Longer MP (multi purpose) ramps have 4 panels. MP ramps also are available with an Edge Barrier, a 2 3/4" tall side panel to eliminate rolling off the edge. DecPac Ramps have a maximum load capacity of 660 lbs.


CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

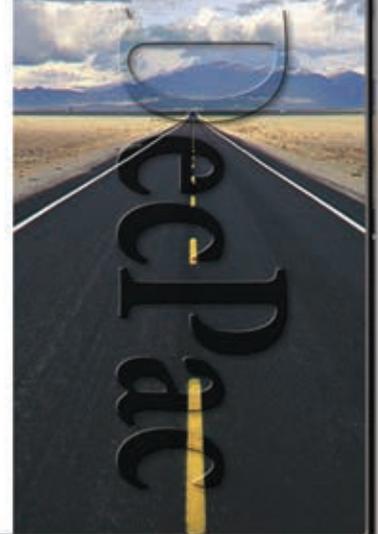

The DecPac Portable Ramp

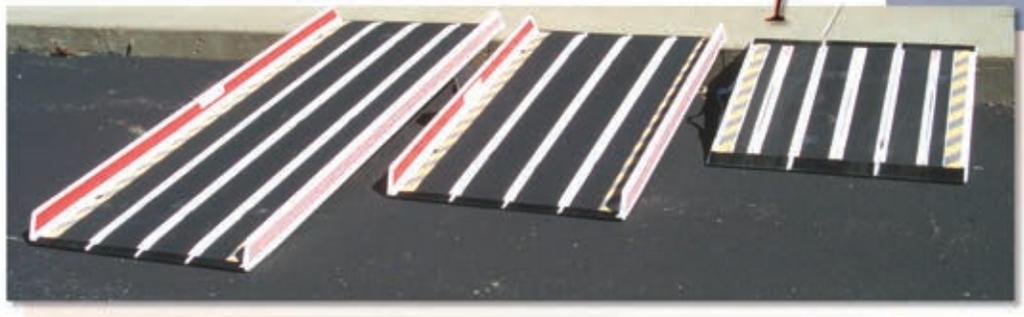

FOLDED SPECIFICATIONS	Length	Width	Weight	Suggested Rise
MP24	2' 4"	29"	8.5 lbs	Up to 7"
MP210	2' 10"	33"	13 lbs	Up to 9"
MP44	4' 4"	29.5"	22 lbs	Up to 13"
MP55	5' 5"	29.5"	26.5 lbs	Up to 16"
MP66	6' 6"	29.5"	31 lbs	Up to 20"
MPEB44	4' 4"	32"	26.5 lbs	Up to 13"
MPEB55	5' 5"	32"	31 lbs	Up to 16"
MPEB66	6' 6"	32"	37.5 lbs	Up to 20"
Capacity	660 lb			

Specifications

To determine ramp length, measure the rise, the distance from the top step to ground. Note, the longer the ramp, the less grade to overcome.

Warranty

DecPac ramps have a 1 year limited warranty against manufacturer defects.



CLARKE HEALTH CARE PRODUCTS, INC.

7830 Steubenville Pike • Oakdale, PA 15071

Phone (888) 347-4537 • Fax (724) 695-2922 • www.clarkehealthcare.com

DecPac

BINNI-Sofa™



Participating Senses:



BINNI-Sofa is a wearable beanbag that provides a child with his own personal space. Curl up, relax or wear a **Binni** and have fun rolling around. The multiple openings allow him to step into it in any manner as he sees fit. The **BINNI** provides soft pressure and becomes a 'soft belly' just calling for a warm hug from mom and dad.

The outer, stretchable, polyester fabric cover is removable and washable. An inner bag contains the "stuffing" beads. Rated for ages 3 and up.

HxWxL		12" x 19" x 17"
Weight		3 lbs.
Part Number	Green	MS10040001GR
	Orange	MS10040001OR

Warranted for 1 year against defects in material or workmanship, from date of purchase. Excludes reasonable wear and tear.

Participating Senses:



MISHU Basic is a weighted doll that provides deep pressure stimulation. Designed to fit shoulders and neck, it can be used in the classroom, during meals or doing your homework. Deep pressure promotes focus after activity.

Mishu's weight can be adapted to the child's age, wishes and needs by inserting or removing the weight pads, 26.5 oz and 9.5oz, in the zippered pockets. The soft polyester fabric can be laundered after removing weights.

HxW	16" x 12"
Weight	5 lbs.
Part Number	MS20010001

Warranted for 1 year against defects in material or workmanship, from date of purchase. Excludes reasonable wear and tear.



SafetyGrip® Bedside Balance Assist Rails

Rails and bars used for assistance while moving, repositioning or sitting up in bed. Assist rails provide stability to balance when standing or transferring from bed.



HER100 SafetyGrip® is a 360° rotating assist rail with both vertical and horizontal grab areas to assist transfer onto and off the bed. Rotation locks every 90°. Using various sections, rail is height adjustable from 34.3"-44.12", can be safely used on all mattress thicknesses and sizes. Vertical and horizontal grab areas with SafetyGrip® coating are 12" long.

Made of steel with white, powder-coat finish, that is anti-microbial and anti-bacterial. Mounts to beds with L-shape angle iron frame. 300 lb. capacity, weighs 19 lbs. 1 year warranty



HER450 SafetyGrip® Balance-Assist bedside rail with multiple height adjustment options. Straps around box spring or frame, fits twin to king size beds, mattress height standard to pillow-top. Top of bar sits 14 ¾" above bed frame, additional cross bar adds 5.25" height, standing 20" above the bed frame or box spring. Grab area is 18.25" with patented SafetyGrip® surface.

Made of steel with white, powder-coat finish, that is anti-microbial and anti-bacterial. 300 lb. capacity, weighs 8.3 lb. 1 year warranty



HSG450 SafetyGlo® balance-assist bedside rail with patented SafetyGrip® grab surface which gives a great grip every time. The unique gentle glow assists the user to locate the grip, gain balance and gather their bearings. Glow is energized by sunlight or room light. Grab area is 18 ¼". Rail attaches with adjustable traps around bed frame or box spring. Top of bar sits 14 ¾" above bed frame, additional cross bar, included in package adds 5.25" height, standing 20" above the bed frame or box spring.

Made of steel with white, powder-coat finish, that is anti-microbial and anti-bacterial. 300 lb. capacity, weighs 8.3 lb. 1 year warranty

These products are not intended to be used as bed restraints to keep someone from falling out of bed.




70002 M-Rail® is an easy bedside mobility solution, fitting snugly against your mattress and helps you get comfortable, reposition, sit-up, balance, stand or transfer to a chair or wheelchair. M-Rail® fits all mattress sizes. Lightweight, compact, and ideal for both home and travel. Adjustable height from 14.5" to 20" above bed frame. Compact and lightweight, its contoured padded handrail provides a soft, secure grip.

Easy to assemble with no tools required. Made of steel, M-Rail® has a 300lb. capacity, weighs 7 lbs. 1 year warranty.



70000 Arcorail® is a rotating, professional-grade bedside handrail, which provides maximum bedside mobility for home, clinical care, and hospital settings. Patented rail rotates 360° and locks in place every 90°; totaling four locked positions. Vertical floor tube has multiple height adjustments, 23"-33" and fits all bed sizes. Second loop can be removed for taller mattresses up to 24" high. Arcorail®'s surface is treated with an anti-bacterial powder coating.

All components are included to attach securely to a variety of bed frame styles, L-shape, tubular, square or rectangular channels. Arcorail® has a 300 lb. capacity, weighs 20 lbs. 1 year warranty.

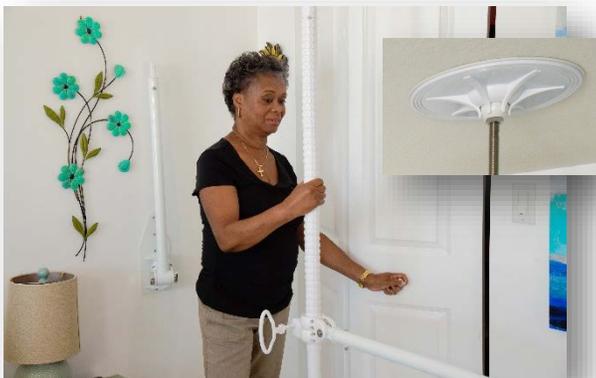

HER611 MoveEasy® Pole is a free standing vertical pole that securely fixes from the floor to the ceiling, anywhere in a room, as long as it has a regular drywall ceiling attached to wooden or steel studs. Floor plate can be screwed into floor for heavy, aggressive use. Easy to store and put away, its fast removal feature in high-occupancy and rotating facilities makes it ideal for facility or home.

Assembles in sections, steel pole adjusts from 79.5" – 109.5" with 39.25" of SafetyGrip®, patented grab surface. Unlike competitors poles which ship as one 8' box, longest section is less than 36" for easy transport or storage. 200 lb. capacity, 1 year warranty.

Photo shows optional horizontal bars not included in HER611.

Model	Width of grab area	Height above mounting point	Capacity	Mounting style
HER100	12"	34.3"-44.12"	300 lb.	Mounts to L-shaped angle iron bed frame
HER450	18.25"	14.75"/20"	300 lb.	Strap around box spring or frame, twin to king
HSG450	18.25"	14.75"/20"	300 lb.	Strap around box spring or frame, twin to king
70002	18"	14.5"/20"	300 lb.	Strap around box spring or frame, twin to king
70000	14"	23"-33"	300 lb.	Mounts to bed frame; L-shaped angle iron, round tube, square or rectangular channel frames
HER611	39.25"	79.5"-109.5"	200 lb.	Pressure fits floor to ceiling

These products are not intended to be used as bed restraints to keep someone from falling out of bed.